

PHILIPS

Heated Humidifier

DreamStation Go

Thinnest travel humidifier¹

Uses tap or bottled water³

Water-saving technology

Works with any mask you like²

HH1531/05



Smart humidification. All night long.

Don't give up the comfort of heated humidification when you're on the go. With its water-saving technology, the convenience of using tap water, and a compact design, you can enjoy heated humidification anywhere your travels take you.

Heated Humidifier

- Use any mask you like
- Freedom to use tap water³
- Manages water usage evenly for comfort throughout the night
- Smart humidification. All night long.
- Lightweight, compact design makes it easy to pack and carry


Highlights

Works with any mask



When you've found a mask you're comfortable with, the last thing you want to do is be forced to switch your mask while travelling. The DreamStation Go heated humidifier is designed to work with any mask that works for you, so you can enjoy on-the-go heated humidification from the comfort of your own mask.²

Freedom to use tap water³

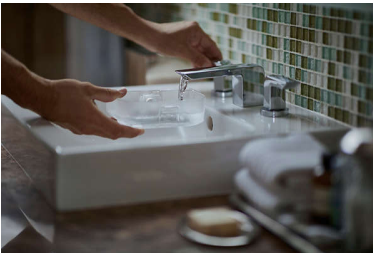

With the freedom to use tap, bottled or distilled water³, the DreamStation Go heated humidifier reduces the inconvenience of taking along heavy, bulky containers or searching for distilled water when you arrive. That means you can choose to travel light and still enjoy the same comfort of heated humidification that you experience at home.

Water saving technology



The DreamStation Go heated humidifier's water-saving technology adapts to a wide range of ambient conditions to deliver humidification evenly throughout the full night's sleep.⁴ You can enjoy the same all-night comfort you experience at home, no matter where your travels take you.

Monitors and adjusts therapy

DreamStation Go features water-saving technology designed to give you a full night of heated humidification⁴ by actively monitoring ambient conditions and therapy settings and automatically adjusting moisture delivery when humidifier target time is enabled. Also, In the event you underfill or forget to fill the humidifier, our heater plate shutdown feature will take over to prevent dry, uncomfortable air from being delivered to you for the remainder of the night.

Lightweight design

The DreamStation Go heated humidifier is designed to make packing a breeze. It's the thinnest travel heated humidifier among leading PAP brands⁵, so it's easy to take along whether your traveling by land, sea, or air.

Specifications

Humidifier Details

Size: 14.7cm x 5.9 cm x 15.1 cm

Weight: 344.7 g

Humidifier Capacity: 240 ml

Tubing compatibility: Micro-flexible 12 mm, 15 mm, and 22 mm standard tubing

Water compatibility: Distilled, tap, bottled*

Humidifier usage while on Battery power: No

Warranty: 2 years



* 1. When compared against Transcend, HDM Z1, AirMini (does not offer heated humidifier), Devilbiss
 * 2. DreamStation Go is equipped with a universal adapter designed to fit all standard 22mm OD connections.
 * 3. In very hard water areas, it is recommended that you use distilled or bottled water.
 * 4. Adjusts humidification output when water-saver target time is enabled. Refer to the User Manual for more information.
 * 5. When compared against ResMed AirSense, AirStart; Devilbiss IntelliPAP 1 and IntelliPAP 2; Fisher and Paykel – SleepStyle, Icon