

SpinVision

Rotational IVUS pullback device

Accurate. Reliable. **Cost effective**.

Accurate

- Long 15 cm, smooth automatic pullbacks facilitate accurate vessel and lesion measurements
- LED display provides pullback position information with accuracy to 0.1 mm
- Maximum pullback accuracy facilitates clinical research and trials

Reliable

- Single board design streamlines device communication to the system
- 30° spring lock design confirms solid catheter to pullback device connection
- Robust pigtail power connector cable design facilitates consistent pullback device to system communication


Cost effective

 Re-usable sled design eliminates need to purchase, stock and track disposables



Length	36.6 cm (14.4")
Width	10.2 cm (4.0")
Height	9.9 cm (3.9")
Weight	1.3 kg (2.9 lb) or 1.5 kg (3.2 lb)
Pullback length	15 cm
Pullback speed options	0.5 & 1.0 mm/sec



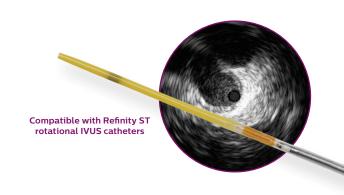

Ordering information

SpinVision option	891 SV
SpinVision sterile drapes (order by each)	891 BG
Refinity ST 45MHz rotational IVUS catheter	Ref 89900 (Includes one sterile drape)

3721 Valley Centre Drive, Suite 500 San Diego, CA 92130 USA www.philips.com/IGTdevices