

# PHILIPS

IPL hair removal  
system

Lumea

For use on body only

15 minutes to treat lower legs

Lifetime >100.000 light  
pulses

SC2001/00



## Prevent hair regrowth

### IPL hair removal: enjoy smooth skin every day

Philips Lumea IPL hair removal system works wonders to prevent the reappearance of hair all over the body. Through gentle pulses of light, applied regularly, Philips Lumea IPL keeps skin silky-smooth every day.

#### **Long lasting smoothness**

- Effortlessly effective
- Enjoy smooth skin every day

#### **Hair removal with IPL technology**

- Adapted for safe use in the comfort of your own home
- Suitable hair and skin types

#### **Gentle and effective, even on sensitive skin**


- Five adjustable light energy settings
- Application areas

#### **Maximum freedom, minimum hassle**

- Cordless convenience
- No replacement parts needed

# Highlights

## Long lasting smoothness

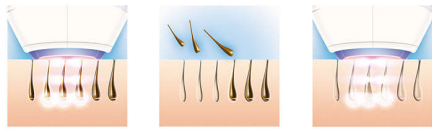

Philips Lumea applies gentle pulses of light to the hair root. As a consequence the hair sheds naturally and hair regrowth is inhibited. Repeating this treatment regularly will keep your skin smooth everyday.

## Effortlessly effective



Our clinical studies have shown significant hair reduction in as little as two bi-weekly treatments, with optimal hair-free results achieved after four to five bi-weekly treatments. To maintain these results, simply repeat treatment every few weeks. Clinical studies show results are best maintained when the treatment is repeated every two to four weeks; however the time between intervals can be extended depending on your hair's regrowth.

## IPL technology



Philips Lumea uses an innovative light-based technology called IPL (Intense Pulsed Light), derived from the technology used in

professional beauty salons. Philips has adapted this technology for safe and effective use in the comfort of your own home. Philips worked closely with leading dermatologists to develop our breakthrough hair removal system. For more than 10 years we carried out extensive consumer research with more than 1,500 volunteers.

## Gentle application



Philips Lumea has five adjustable light energy settings to ensure gentle but effective treatment. Used correctly, Philips Lumea's IPL treatment is safe and gentle to use even on sensitive skin and sensitive body areas.

## Cordless convenience



Lumea's unique cordless design is perfect for treating the areas you want, when you want.

## No replacement parts needed



Philips Lumea is a complete solution that comes ready to work right out of the box. More importantly, it does not require any replacement lamps or gels.

## Application areas



Use with confidence to remove hair on body areas below the neckline: legs, underarms, bikini area, arms and stomach.

## Suitable hair and skin types



Philips Lumea works effectively on (naturally) dark blonde, brown and black hairs. As with other light-based treatments, Philips Lumea is NOT effective on red, light-blonde or white/grey hair. Philips Lumea is also NOT suitable for very dark skin.

# Specifications

## Power

**Voltage:** 100-240 V

## Service

**Warranty:** 2 years global warranty

## Technical specifications

**High performance lamp:** Delivers <gt/>  
50,000 flashes

## Application areas

**Body areas:** Legs, Arms, Armpits, Belly, Bikini

## Attachments

**Body attachment (3cm2):** For use below the neckline

## Items included

**Charger:** Multi voltage charger

**Instructions for use:** Demo DVD, Quick Start Guide, User manual

**Storage:** Pouch

## Power

**Charging time:** 90 min

**Running time:** Min. 150 flashes

## Safety and adjustable settings

**Integrated UV filter:** Protects skin from UV light

**5 light energy settings:** Adjustable to your skin type

**Integrated safety system:** Prevent unintentional flashing

## Application time

**Lower leg:** 8 min

**Armpit:** 1 min

**Bikini line:** 1 min

## Application mode

**Step and Flash:** For treatment on curvy areas

**Corded / cordless use:** Cordless

