

Philips Lumea IPL 7000
Series
IPL Hair removal device

5 manual intensity settings

2 attachments: body, face
Lumea IPL App
Corded use

BRI921/60



Great results, great value

Enjoy 12 months of hair-free smooth skin*

Get hair-free smooth skin with our Lumea IPL 7000 Series. Treatments are gentle and effective with attachments for each body area.

Gentle and effective for long-lasting smooth skin

- Treat only 2x a month for fast results
- Developed with dermatologists to be easy and effective
- Extra-long cable for added flexibility while treating

Full solution for face and body

- Full solution for face and body with 2 attachments
- Comes with Satin Compact pen trimmer
- Suitable for a wide range of skin tones and hair colors

Personalized IPL hair removal treatment

- Skin Tone Sensor and 5 intensity settings
- Optimize your routine with the Philips Lumea IPL app

PHILIPS

IPL Hair removal device
5 manual intensity settings 2 attachments: body, face, Lumea IPL App, Corded use

BRI921/60

Highlights

2x a month for fast results



Every 2 weeks to get started — that's half as many treatments as other brands. Followed by touch-ups just once a month. That's it. Covering both lower legs takes 15 minutes.

Skin Tone Sensor



Choose from 5 intensity settings for a comfortable experience. The Skin Tone Sensor prevents you from treating areas of your skin that are too dark for IPL treatment.

Extra-long cable




Convenient to use thanks to an extra long cable for easy access and enhanced maneuverability.

2 dedicated attachments




Different heads for body and face work with speed and precision for gentle use, even on sensitive areas. The face attachment has an integrated light filter for precise treatment.

Developed with dermatologists

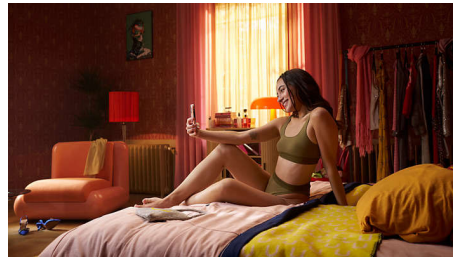

As a leader in health technology, Philips developed Lumea IPL in consultation with dermatologists, for easy and effective use in the safety of your home. Derived from technology used in professional salons, Lumea IPL offers gentle treatment, even on sensitive areas.

Philips Lumea IPL app

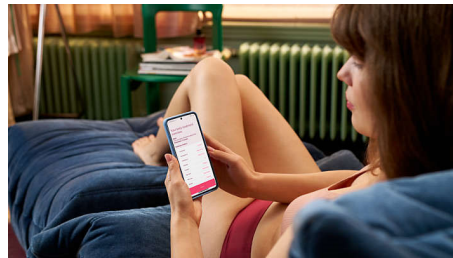

Our free coaching app helps you plan and stick to your treatment schedule, then takes you

through each session step by step. Downloaded by more than 2.1 million users.

Satin Compact pen trimmer



Complete your beauty routine: The Satin Compact pen trimmer is a discrete on-the-go beauty tool to quickly and easily remove even the finest facial hairs.

For many skin and hair types



IPL needs contrast between the pigment in the hair color and the pigment in the skin tone, therefore works on naturally dark blonde, brown and black hair and on skin tones from fair to medium brown (I-IV).

Specifications

Safety and adjustable settings

- Integrated UV filter: Protects skin from UV light
- 5 light energy settings: Adjustable to your skin type
- Integrated safety system: Prevent unintentional flashing

Application mode

- Slide and Flash: For quick application
- Corded / cordless use: Corded use
- Stamp and Flash: For treatment on small areas

Pen trimmer

- Accessories: Eyebrow comb, Cleaning brush
- Ease of use: Compact handle
- Power: AAA battery

Technical specifications

- High performance lamp: Built to last, 250,000 flashes, equivalent to 20 years of lamp lifetime**
- Voltage: 100-240 V

Service

- Warranty: 2 years global warranty + 1 extra year of warranty upon product registration within 90 days

Items included

- Instructions for use: User manual
- Storage: Beauty pouch
- Adapter: 24V / 1500mA

Application time

- Lower leg: 15 min
- Underarms: 2,5 min
- Bikini line: 4 min
- Face areas: 2 min

Technical specs. attachments

- Body attachment: Window size: 4.0 cm², Treatment for body: legs, arms, stomach and underarms
- Face attachment: Window size: 2 cm², extra filter, Treatment for face: upper lip, chin and jawline



Issue date 2024-03-24

Version: 8.8.1

© 2024 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* Median hair reductions after 12 treatments: 82% on lower legs
** ** When following treatment schedule. Calculated for use on lower legs, bikini, armpits and face. The lamp lifetime does not extend the Philips 2 years worldwide guarantee