

PHILIPS

Infrared lamp

300 W

For deep, zone treatment

Adjustable angle



PR3120/00



Effective*, drug-free pain relief

Clinically proven technology*

The InfraCare technology provides comforting warmth by penetrating deep into the skin. This increases local blood circulation & enhances the supply of oxygen & nutrients to help treat sore muscles and effectively recover from joint pain

Effective* pain relief

- Pain relief and better mobility
- Helps the body to heal in a natural manner

Proven health benefits

- Clinically proven
- Classified medical device

Zone treatment

- Localized pain relief

Highlights

Effective* pain relief




Our clinically proven technology provides relief from muscle and joint pains without the use of drugs. The InfraCare PR3120 stimulates local blood circulation and increases the local supply of oxygen for an effective zone treatment. This in turn provides a relaxing affect to the injured area and helps in an effective* recovery from pain.

Deep infrared treatment



InfraCare is designed to deliver the highest level of infrared waves reaching the deepest levels of the skin. Infrared warmth increases the local blood circulation, which speeds up the removal of metabolic end products and carbon dioxide from muscle tissue.

Clinically proven

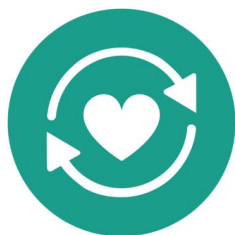

Trusted technology used by professionals. Clinically proven to objectively lead to functional joint improvement*. InfraCare has proven to lead to significant pain reduction in 71% of rheumatoid arthritis patients and in 72% of low back pain patients*.

Classified medical device



A safe and effective way to treat your pain. InfraCare complies with IEC 60601-1 3rd edition for medical devices.

Adjustable angle



The InfraCare 300W Lamp can be easily positioned from 0-40 degrees to relieve pain in a comfortable and an effective* manner.

Specifications

Easy to use

Adjustability: 0-40 degrees (backwards)
1 handle + 1 handgrip
Mechanical timer

Focused treatment

Treatment area: 40 x 30 cm
Infrared halogen lamp: 300 W

Logistic data

Country of origin: Hungary

Medical appliance

Medical Device Directive: 2007/47/EC, MDD 93/42/EEC

Safety

IEC certified: IEC 60601-1, 3rd edition
30 minutes auto shut-off

Technical specifications

Power: 300 W
Voltage: 110 V or 220-230 V
Frequency: 50-60 Hz
Insulation: Class II (double isolation)
Lifetime of lamp: 500 hour(s)
Cord length: 250 cm

Weight and dimensions

Product dimensions: 18,0 x 20,6 x 31,6 cm
Product weight: 2 kg
F-box dimensions: 22,3 x 24,9 x 33,4 cm

