

PHILIPS

Infrared muscular pain relief

200W

Focused treatment area



HP3621/02



Effective pain relief

Deep penetrating infrared warmth

Philips InfraCare gives effective pain relief for muscles and joints. Its 200 Watt infrared warmth is very comfortable and penetrates deep into skin, improving blood circulation thereby giving effective relief for focused areas of 20x30 cm.

Effective pain relief

- Deep penetrating Infrared warmth

Extra focused treatment

- 200 Watt infrared halogen lamp

Easy to use


- Adjustable angle



asimpleswitch.com

Highlights

Deep penetrating Infrared warmth

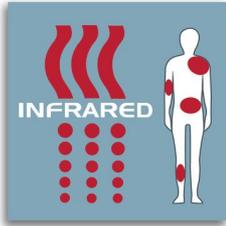

Infrared light has been proven to effectively relieve aches and pains caused by muscular pains and stiff joints. The comfortable warmth created by infrared light penetrates deep into your skin, stimulates blood circulation and warms your muscles. As your muscles are soothed, they automatically loosen up and relax. Because it makes the tissue more pliable, warmth also reduces stiffness and makes joints more flexible.

200 Watt infrared halogen lamp



The Philips InfraCare HP 3621 has been optimised for effective relief of pains in areas covering 20 x 30 cm, for example shoulder,

elbow, calf, or neck. The InfraCare features innovative infrared halogen lamp technology. The special optics, filter and the powerful halogen lamp of 200 Watt have been optimised to provide a clear focused treatment as the infrared warmth is more equally divided over the treatment area.

Adjustable angle



The infrared lamp can easily be positioned upto 40 degrees backwards.

Specifications

Technical specifications

- Power:** 200 W
- Voltage:** 110 or 220 or 230 V
- Frequency:** 50/60 Hz
- Lifetime of lamps:** 500 hour(s)
- Cord length:** 2 cm
- Insulation:** Class II (double isolation)

Weight and dimensions

- Product dimensions:** 20x30.9x17.2 (WxHxD) cm
- Product weight:** 1.4 kg
- F-box dimensions:** 28x38.6x24.5 (WxHxD) cm
- F-box weight:** 1.8 kg
- No. of F-boxes in A-box:** 3 pcs
- A-box dimensions:** 29x40.4x74.5 (WxHxD) cm
- A-box weight:** 5.4 kg
- Qnt. on Euro pallet:** 60 pcs

Logistic data

- CTV code:** 884362102000
- Country of origin:** Germany

Safety

- IEC certified:** In line with IEC 60601 and 60335

Easy to use

- Adjustability:** 0-40 degrees (backwards)
- Soft-touch handgrip:** for easy transport and positioning
- On/off switch**

Medical appliance

- Proven health benefits:** supported by Medical Device Directive certificate

Effective pain relief

- Deep penetrating warmth**
- Even warmth distribution:** Now, upto 3 x more comfortable warmth

Focused treatment

- Treatment area:** 20x30 cm
- Infrared halogen lamp:** 200 W



asimpleswitch.com

Philips Green Logo

Philips Green Products can reduce costs, energy consumption and CO2 emissions. How? They offer a significant environmental improvement in one or more of the Philips Green Focal Areas – Energy efficiency, Packaging, Hazardous substances, Weight, Recycling and disposal and Lifetime reliability.

