

# PHILIPS

IPL Hair removal device with SenseIQ

Lumea IPL 9000 Series

SmartSkin Sensor

4 intelligent attachments: body, face, bikini, underarms

Lumea IPL App

Both cordless and corded use



BRI957/60



## Cordless ease, personalized treatment with SenseIQ

Enjoy 12 months of hair-free smooth skin\*

Enjoy hassle-free, personalized treatment with the cordless version of our fastest Lumea IPL device. SenseIQ technology with smart attachments and guidance in the Lumea IPL app to give you long-lasting smooth skin.

### Specialized treatment of different body areas

- Optimize your routine with the Philips Lumea IPL app

### Gentle and effective for long-lasting smooth skin

- Treat only 2x a month for fast results
- Developed with dermatologists to be easy and effective
- The only Lumea with cordless option

### Full solution for face and body

- Full solution for face and body with 4 smart attachments
- Suitable for most skin tones and hair colors

### Personalized IPL hair removal treatment

- Gentle and comfortable with SenseIQ


# Highlights

## 4 smart attachments




Specially designed attachments perfectly fit the curves of your body and trigger the most effective programs for each body area when they are connected. Face: Flat design and small window with UV filter. Body: Curved inward with a large window. Bikini and Armpit: Curved outwards for hard-to-reach areas.

## Philips Lumea IPL app

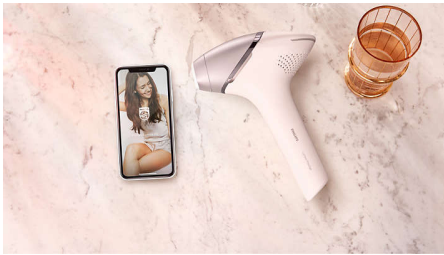

Our free coaching app helps you plan and stick to your treatment schedule, then takes you through each session step by step. Downloaded by more than 2.1 million users.

## Developed with dermatologists

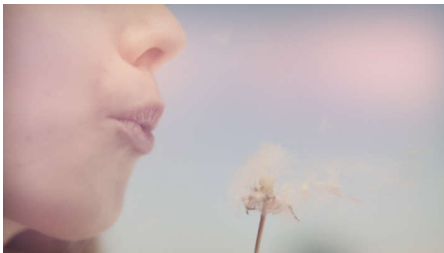

As a leader in health technology, Philips developed Lumea IPL in consultation with dermatologists, for easy and effective use in

the safety of your home. Derived from technology used in professional salons, Lumea IPL offers gentle treatment, even on sensitive areas.

## 2x a month for fast results



Every 2 weeks to get started — that's half as many treatments as other brands. Followed by touch-ups just once a month. That's it. Covering both lower legs takes 8.5 minutes.

## SenseIQ technology



Lumea 9000 Series has five easily adjustable light settings. Our SmartSkin sensor reads your skin tone and helps you find the most comfortable setting. Intelligent attachments adapt treatment for each body area.

## Cordless if you want



For extra convenience, Lumea 9000 Series is our only IPL with cordless functionality, which makes it even easier to reach anywhere with precision and to use it from anywhere you want.

## For most skin and hair types



IPL needs contrast between the pigment in the hair color and the pigment in the skin tone, therefore works on naturally dark blonde, brown and black hair and on skin tones from fair to dark brown (I-V).

# Specifications

## Technical specs. attachments

**Body treatment attachment:** Shape: Convex curved, Window size: 4.1 cm<sup>2</sup>, Tailored treatment for body: legs, arms and stomach

**Face treatment attachment:** Shape: Flat, Window size: 2 cm<sup>2</sup>, extra filter, Tailored treatment for face: upper lip, chin and jawline

**Bikini treatment attachment:** Shape: Concave curved, Window size: 3 cm<sup>2</sup>, extra filter, Tailored treatment for sensitive bikini area

**Underarm treatment attachment:** Shape: Concave curved, Window size: 3 cm<sup>2</sup>, Tailored treatment for underarms

## Safety and adjustable settings

**Integrated UV filter:** Protects skin from UV light

**5 light energy settings:** Adjustable to your skin type

**Integrated safety system:** Prevent unintentional flashing

**SmartSkin sensor:** The right setting on demand

## Application mode

**Slide and Flash:** For quick application

**Corded / cordless use:** Both cordless and corded use

**Stamp and Flash:** For treatment on small areas

## Technical specifications

**Battery type:** Rechargeable battery

**Voltage:** 100–240 V

**High performance lamp:** Built to last, 450,000 flashes, equivalent to 39 years of lamp lifetime\*\*

## Service

**Warranty:** 2 years global warranty + 1 extra year of warranty upon product registration within 90 days

## Items included

**Instructions for use:** User manual

**Storage:** Luxurious pouch

**Adapter:** 19.5V / 4000mA

## Application time

**Underarms:** 2,5 min

**Bikini line:** 2 min

**Face areas:** 1,5 min

**Lower legs:** 8,5 min

