

PHILIPS

Compressor
Nebulizer System

InnoSpire Essence

Compact home-use device

Fast, effective treatment

Nebulizer kit included



HH1340



A nebulizer system you can rely upon

When reliability is key choose InnoSpire Essence. Essence provides fast, efficient nebulized drug therapy for you and your family at home.

Compact design

- Unique diamond jet design

Fast treatment times

- Boosts air flow for a faster treatment time

Proven SideStream technology

- Proven nebulizer technology
- Robust design meets reliability

Short treatment times

- Treatment delivered in 6-8 minutes*
- Treatment delivered in 6-8 minutes*

For use in the home

- Designed for intermittent use
- Compatible with commonly prescribed nebulizer medications


Highlights

Treatment time



InnoSpire Elegance combined with the SideStream Nebulizer, is designed to provide fast drug delivery and short treatment times.

For home use

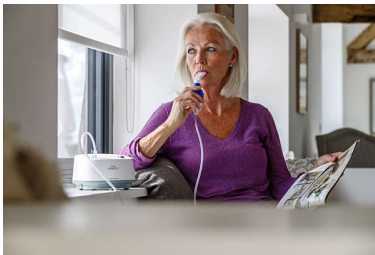

Designed for intermittent use and suitable for the home, InnoSpire Essence provides the consistent delivery high-quality aerosol treatment that you can rely on.

SideStream nebulizer system



Together with the SideStream reusable nebulizer, InnoSpire Elegance provides consistently fast and efficient aerosol therapy.

Faster drug delivery

InnoSpire Essence is a compact and reliable compressor nebulizer system suitable for frequently nebulizing most commonly prescribed aerosol medications.

Versatility



Capable of nebulizing a wide variety of medications that are routinely with nebulizers used to treat various respiratory conditions.

Treatment time

InnoSpire Essence, combined with our SideStream nebulizer, is designed to provide fast drug delivery and short treatment times.

Proven nebulizer technology



Reduced Treatment Times - Our Sidestream nebulizer is specially designed to keep your treatment times as short as possible. As the InnoSpire Essence pushes air through the tube, additional air is drawn in through the top of the nebulizer to create a fast flow of aerosol medication for you to inhale. Treatments can take as little as 6 minutes* Consistent delivery of high quality aerosol therapy - Sidestream is designed to withstand frequent and repeated use and will continue to efficiently deliver your nebulizer medication until it is time for replacement.

Power and style all in one



InnoSpire Essence is designed to both fit your lifestyle and cope with frequent use. You can be confident that you will receive high quality aerosol for the life of the nebulizer.

Specifications

Service

Warranty: 3 years

Design

Compact: Weighs only 1.5 kg

Operating parameters

Electrical rating/ power: 230 V / 50 Hz, 1 Amp

Maximum flow rate: 8 LPM

Mode of operation: Intermittent (30 min. on/off)

Average flow rate: 6 LPM @ 10 psi

Product details

Weight: 1.5 kg

Size: 165x165x108 mm

Noise level: 58 ±3 dBA

Treatment time (6-8 minutes*): * 2.5ml Salbutamol dose

Performance data

MMAD: 2.9 µm

Fine particle fraction: 77% below 5 µm

Nebulization capacity: 8 ml

