

# PHILIPS

## Auto CPAP

### DreamStation 2 Advanced

Includes Ramp Plus

Color touchscreen

Integrated humidifier

3 adjustable comfort settings

HH1551/01



## Advanced care designed to meet your needs

### Includes 12mm Heated tubing

DreamStation 2 Auto CPAP Advanced is designed with your comfort and ease in mind. Now with Ramp Plus, you can pick a starting pressure that is comfortable for falling asleep. It also features a colorful touchscreen for easy navigation.

#### **DreamStation 2**

- Designed for comfort
- Easy to use
- Quiet sound
- Small size
- Keeps you motivated

# Highlights

## Comfort




With Ramp Plus, you have the ability to pick the starting pressure that is comfortable to fall asleep to and will adjust should an event occur.

## Easy



Three comfort settings are designed to be easy to adjust while therapy is running. Integrated connectivity within the device allows you to stay connected with your care team. 93% of users describe the CPAP as “easy to use”.<sup>1</sup>

## Quiet

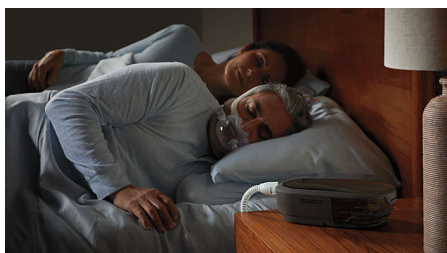

Users rated DreamStation 2 Auto CPAP Advanced quieter than their prescribed device.<sup>2</sup>

## Small - 12mm tube



DreamStation 2 Auto CPAP Advanced is 22% smaller and 16% lighter than the other leading CPAP brands.<sup>3</sup> The device also features a small, micro-flexible 12 mm tube that weighs less than the leading brands' 15mm or larger CPAP heated tubing.\*

## DreamMapper



DreamMapper mobile app allows you to track your progress, set goals and get helpful tips. People who struggled with sleep apnea therapy had a 283% higher success rate when they used DreamMapper.<sup>4</sup>

# Specifications

## Product details

**Weight:** 2.3 lbs without battery; 2.94 lbs with battery

**Dimensions:** 10.11 in L x 6.24 in W x 3.34 in H

**Package contents:** Auto CPAP device, 1 Disposable Ultra-fine filter, 1 Reusable pollen filter, 12mm Tubing, Carry bag, Power cord, Power supply, User Manual

**Tubing compatibility:** 12mm standard/heated, 15mm standard/heated, 22mm standard

## Technical Specifications

**Pressure range:** 4 to 20cm H<sub>2</sub>O

**Flex Pressure Relief:** Off, 1, 2, 3

**Ramp time:** 15 minutes, 30 minutes, 45 minutes

**Humidification:** Adaptive humidification and heated tube settings (optional)

**Altitude compensation:** Automatic

**Filters:** Reusable pollen; disposable ultra-fine

**Device controls:** LCD color touchscreen; push therapy button

**Electrical requirements:** 100-240 VAC, 50/60 Hz, 1.0-2.0 A

**Modes:** CPAP, Auto CPAP



\* 1. Based on 2020 patient preference study with N=444

\* 2. Based on users' rating of sound and comparison from 30-day 2020 patient preference study with N=444. Prescribed devices included Philips DreamStation and ResMed AirSense 10.

\* 3. Based on 2020 study (n=160 patients) where patients tested DreamStation 2 devices that have either prescribed Philips DreamStation devices or ResMed AirSense10 devices.

\* 4. 46% of DreamMapper patients went on to achieve adequate adherence within 90 days versus 12% for the Standard Care group, a relative increase of 283.33%. In a retrospective review conducted by Philips Respironics of the EncoreAnywhere database (white paper) of struggling 19 patients (n=24,176).

\* \*Leading brands: ResMed ClimateLine Air 15mm heated tube, Fisher & Paykel ThermoSmart 19mm heated tube

\* 6. Leading brands: ResMed ClimateLine Air 15mm heated tube, Fisher & Paykel ThermoSmart 19mm heated tube.