

PHILIPS

EnergyLight



HF3308



Feel as energetic

as on a sunny day

Special light device that gives you the daylight—and energy—you need, whatever the season, in a way that easily fits into your daily routine.

Feel energetic in 30 minutes a day

- Original Philips EnergyLight lamps

Easy positioning for personal preference

- Adjustable angle

Safe to use

- UV-free

Proven health benefits


- Supported by the Medical Device Directive (MDD)



asimpleswitch.com

Highlights

Original Philips EnergyLight lamps

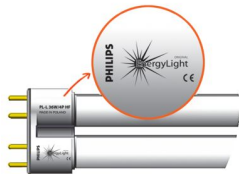

The EnergyLight features original Philips EnergyLight lamps producing light with all the intensity and quality of natural daylight. In fact, it can give you up to 10,000 lux, but without any unpleasant dazzle or glare. If you use the EnergyLight for just 30 minutes a day (there's no need to stare into the lamp all the time) at a light intensity of 10,000 lux, you'll notice an improvement in your energy level and mood after just 5 days' use.

Adjustable angle



The appliance can be positioned at various angles for comfortable use.

UV-free



The EnergyLight uses the full natural daylight spectrum. It meets European safety standards and features a UV-filtering screen. The EnergyLight is therefore safe to use.

Supported by the Medical Device Directive (MDD)



This device has proven health benefits, as certified by the Medical Device Directive (MDD).

Technical specifications

Power: 72 W
Voltage: 220/230
Frequency: 50/60 Hz
Type of lamps: 2 PL-L 36W Philips EnergyLight
Lifetime of lamps: 10.000
Type of ballast: HF-K 236TLD (High Frequency)
Ballast operating frequency: 42.000 kHz
Insulation: Class II (double isolation)
Cord length: 3 m

Weight and dimensions

Product dimensions: 33x57x18 (WxHxD) cm
Product weight: 3.5 kg
F-box dimensions (= A-box): 36x62.2x21.5 (WxHxD) cm
F-box weight (= A-box): 4.4 kg
Qnt. on Euro pallet: 30 (UK 39) pcs

Logistic data

CTV code: 884330801000
Country of origin: Germany

Safety

UV-free: No UV comes through the filtering screen
No flickering: due to high frequency ballast avoiding headaches
Cenelec certified

Easy to use

On product user protocol: Indicates how long to use it at a certain distance
On/off button

Feel energetic in 30 min a day

Original Philips EnergyLight: lamps: 10.000 lux at 20 cm, 2.500 Lux at 60 cm

Easy positioning

Adjustable angle: 20 degrees
Compact size: 18x33 (DxW) cm

Satisfaction indicator

94% satisfied users

Medical appliance

Proven health benefits

