

# PHILIPS

## EnergyLight



HF3319



## Feel as energetic

Boosts your energy as on a sunny day

The Philips Energy Light is a special light therapy device that gives you the daylight and energy for you to retain the summer feeling all year long in a way that easily fits your daily routine.

### **Feel energetic in 30 minutes a day**

- 10,000 lux light intensity
- Clinically proven

### **Comfortable and easy-to-use light**

- Light intensity dimmer
- Adjustable angle
- Treatment time indicator

### **Safe to use**

- UV-free



# Specifications

## Wellbeing from light

Fight the winter blues  
Increase alertness  
Increase energy level

## Comfortable light

10,000 Lux  
Even screen illumination

## Easy to use

Dimmer  
Treatment time indicator  
On/off button

## Easy positioning

Adjustable angle: 20 degrees

## Safety

Complies to IEC 60601-1 Ed. 3  
UV-free: No UV or near UV radiation

## Technical specifications

Power: 75 W  
Frequency: 50 Hz  
Type of lamps: 2 PL-L 36 W Philips EnergyLight  
Insulation: Class II (double insulation)  
Voltage: 230 V  
Lifetime of lamps: 1700 hours (10,000 lux)  
Ballast operating frequency: 45 kHz...125 kHz  
Cord length: 200 cm

## Weight and dimensions

F-box dimensions (= A-box): 36 x 62.2 x 21.5 (W x H x D) cm  
F-box weight (= A-box): 4.4 kg

## Logistic data

Country of origin: Germany  
CTV code: 884331901

## Medical appliance

Medical Device Directive: 2007/47/EC, MDD 93/42/EEC

# Highlights

## 10,000 lux light intensity



The EnergyLight produces light with all the intensity and quality of natural daylight. In fact, it can give you up to 10,000 lux, but without any unpleasant dazzle or glare.

## Clinically proven

The Philips EnergyLight is clinically proven to increase your energy level and help fight winter blues.

## Light intensity dimmer




Fight winter blues and get the energy you have in summer while using the appliance at a convenient level of light you need for easy-reading etc. The EnergyLight features a light intensity dimmer so it can easily be adjusted between zero to about 10,000 lux (at 15-20 cm) matching the intensity of your personal preference.

## UV-free

The EnergyLight uses the full natural daylight spectrum. It meets European safety standards and features a UV-filtering screen. The EnergyLight is therefore safe to use.

## Adjustable angle

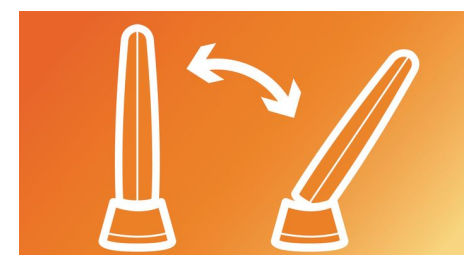

The appliance can be positioned at various angles for comfortable use.

