

# PHILIPS

## energy light

## goLITE BLU

rechargeable

international adapters



HF3330/01



## Feel more energetic

Energy throughout the day, in a natural way

The Philips goLITE BLU is a compact and portable energy light which produces blue light like a clear summer sky. This particular light therapy is known to naturally help increase your energy level whenever you need it and help fight winterblues.

### **Increases energy levels and fights winter blues**

- Uses specific blue light like a clear sunny sky
- Helps to alleviate the symptoms of the winter blues
- Helps to increase energy levels
- Helps to overcome tiredness from jet-lag

### **Comfortable, convenient and easy to use**

- Use for just 15–45 minutes per day to experience results
- Conveniently portable and rechargeable
- Backlit touch screen to adjust timer & intensity
- Extra durable LED lights
- 100% UV free blue light

### **Proven benefits**

- Independent research show improved energy levels with users
- Clinically proven\* to improve mood and energy
- 100 years of Philips expertise in light technology

# Highlights

## Natural blue light



goLITE BLU energy light produces a particular kind of pure blue light that occurs naturally on very sunny, clear days. Special receptors in the eyes absorb this blue light which makes you feel energetic whenever you need it.

## Fights winter blues



In countries with long winters, many people experience loss of energy and low spirits due to lack of light. These symptoms are referred to as the winter blues. They typically start in autumn and may last for several weeks. Philips goLITE BLU is a simple and natural way to alleviate these symptoms and restore your summer mood and energy level.

## Increases energy levels




Our busy lifestyles demand a lot from us and we want to be at our best both during work and private time. Sometimes we experience low moments or dips during the day. Using the goLITE BLU energy light at any time of the day will help you to increase your energy level and alertness in a natural way.

## Helps with jet-lag

Many people who travel regularly experience jet-lag after long flights. Using the goLITE BLU when you feel the lag, gives you a natural energy lift and can help reset your body clock to the new time zone.

## Use for 15-45 minutes per day

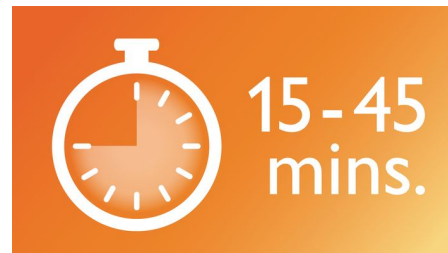

Place goLITE BLU in your field of vision - no need to look directly into the light. Using the goLITE BLU energy light for at least 15 minutes naturally improves your mood and energy level whenever you need it. Using it for longer is absolutely safe. Regular use every morning during winter months will help to alleviate the symptoms of the winter blues.

## Use everywhere



Compact design, storage case, rechargeable battery and international adapters make the goLITE BLU easy to use, easy to store and easy to take with you when you travel.

## Adjustable timer and intensity



Set the timer according to your needs (1-60 min) by using the touch display and choose the right level of brightness out of the four

levels by using the dimmer. Set the treatment time alarm to ensure that you don't miss your session of blue light.

## Long lasting LEDs




The durable LED lights in the goLITE BLU energy light can be used 30 minutes a day for 50 years, without losing its efficacy.

## UV free



The pure, blue light increases your energy naturally, like on a sunny day, but without the UV rays associated with sunlight.

## Improved energy levels

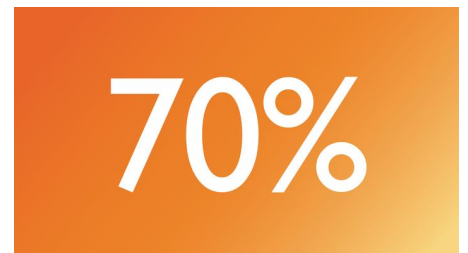

Independent research shows that 70% felt their energy levels improved after 2 weeks use. Bzz Agent UK 2011; n = 823

# Specifications

## Well-being by light

Fight winter blues  
Increase energy level  
Increase alertness  
Fight jet lag

## Comfortable light

Specific blue light  
Light Intensity dimmer: Yes, 4 levels  
Advanced diffusion optics  
Wide treatment field

## Easy to use

Backlit touch screen  
Built in clock  
On/off button  
Pause function  
Treatment reminder alarm  
Treatment timer: 0-60 min

## Portable

Small, compact design  
Internal rechargeable battery  
International power adapters  
Protective pouch included

## Medical appliance

Medical Device Directive: 2007/47/EC, MDD  
93/42/EEC

## Safety and Regulations

CE 0344 Certified  
Complies to IEC 60601-1 Ed. 3  
Not for general illumination  
UV-free: No UV or near UV radiation

## Technical specifications

Power: 12 W  
Voltage: 100-240 V  
Frequency: 50/60 Hz

Insulation: Class II (double insulation)

Type of lamps: LED

Lifetime of lamps: 10,000 hour(s)

Cord length: 180 cm

## Weight and dimensions

Product dimensions: 14x 14 x 2.5 cm

Product weight: 0.4 kg

Box dimensions (WxHxD): 17.8 x 17.8 x 17.8 cm

Box weight: 0.9 kg

Master carton: 2

## Logistic data

Country of origin: China

CTV code: 884333001

