

PHILIPS

sonicare



Sonic electric
toothbrush

ProtectiveClean 4300

Built-in pressure sensor

1 cleaning mode

1 x BrushSync feature



HX6803/03

Cleaner teeth. Gently does it.

Remove up to 7 x more plaque than a manual

Feel the difference of a gentle clean with our Pressure Sensor while removing up to 7 x more plaque than a manual toothbrush

Proven to improve oral health

- Up to 7x more plaque removal than a manual toothbrush
- Safe and gentle on sensitive areas, orthodontics and dental work

Designed around you

- Easy-start programme builds your Philips Sonicare routine
- Encouragement to brush thoroughly
- Makes traveling easier for you
- Always know when to replace your brush heads

Innovative technology

- Lets you know if you're pressing too hard
- Philips Sonicare's advanced sonic technology
- Connects smart brush handle and smart brush heads

ProtectiveClean 4300

Sonic electric toothbrush

HX6803/03

Highlights

Say goodbye to plaque



Densely-packed, high-quality bristles remove up to 7 x more plaque than a manual toothbrush. And the specially curved power tip makes reaching the teeth at the back of your mouth a breeze.

Safe and gentle




You can be sure of a safe brushing experience: our sonic technology is suitable for use with braces, fillings, crowns and veneers, and helps prevent cavities and improve gum health.

The only true sonic technology



Philips Sonicare's advanced sonic technology pulses water between your teeth, and its brush strokes break plaque up and sweep it away for an exceptional daily clean.

Monitor your brushing pressure

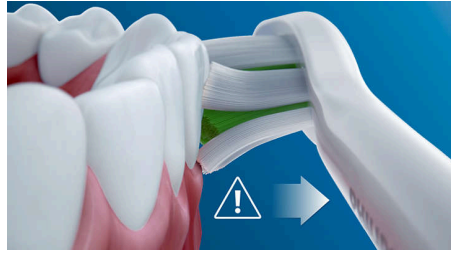

Brushing too hard can damage your teeth and gums. To prevent this, your Philips Sonicare ProtectiveClean emits a gentle pulsing sound to remind you to ease off on the pressure.

BrushSync technology



A microchip-enabled technology that detects and synchronises the smart brush head with the smart handle. The smart handle and smart brush head pair is a powerful combination that enables Smart replacement reminders.

BrushSync replacement reminder




All brush heads wear out over time. But our BrushSync technology tracks how long you've been using your brush head for, and how hard you've been brushing. When it's time to replace it, a light on your handle and a short beep will let you know. That way, you can be sure your brush head is doing a good job.

Ease into your new routine



Starting something new can take some getting used to. Our easy-start programme gives you the option of a gradual, gentle increase in brushing power with the first 14 uses of your new toothbrush.

Helpful timers

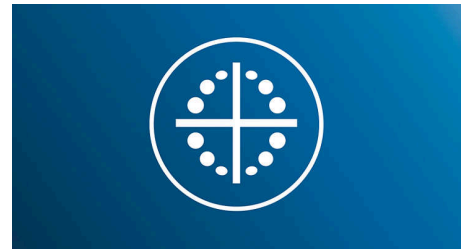

Need an electric toothbrush with a timer? Our QuadPacer lets you know when you've spent just the right amount of time cleaning each part of your mouth, while our Smartimer tells you when you've brushed for the recommended two minutes.

Handy travel case and charger



Our premium travel case lets you store your toothbrush hygienically, while our compact charging base keeps you topped up when you're on the go. You do get two weeks of regular use from a single full charge, but the charger's just the ticket for longer trips.

ProtectiveClean 4300

Sonic electric toothbrush

HX6803/03

Specifications

Power

Voltage: 110-220 V

Technical specifications

Operating time (full to empty): Up to 2 weeks

Battery: Rechargeable

Battery type: Lithium ION

Energy consumption: Standby without display
<0.5 W

Design and finishing

Colour: Light blue

Service

Warranty: 2-year limited warranty

Ease of use

Handle: Slim ergonomic design

Handle compatibility: Easy click-on brush heads

Battery indicator: Light shows battery status

Brushing time: Up to 2 weeks

Items included

Handle: 1 ProtectiveClean

Brush heads: 1 C2 Optimal Plaque Defence

Travel case: 1

Charger: 1

Cleaning performance

Performance: Removes up to 7x more plaque*

Health benefits: Helps reduce cavities

Timer: Quadpacer and Smartimer

Speed: Up to 62,000 brush movements/min

Pressure feedback: Vibrates handle to alert user

Modes

Clean: For exceptional everyday cleaning

2 intensities: Low, High

Smart sensor technology

Pressure sensor: Alerts when brushing too hard

BrushSync Replacement Reminder: Always know when to, replace brush heads

