

Philips
EnergyLight

HF3309/01



Feel as energetic

as on a sunny day

Special light device with intensity dimmer that gives you the daylight—and energy—you need, whatever the season, in a way that easily fits into your daily routine.

Safe to use

- UV-free

Proven health benefits

- Supported by the Medical Device Directive (MDD)

Feel energetic in 30 minutes a day

- Original Philips EnergyLight lamps

Easy positionable for personal preference

- Adjustable angle


Comfortable light

- Light intensity dimmer

PHILIPS

Highlights

Original Philips EnergyLight lamps

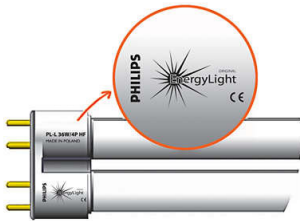

The EnergyLight features original Philips EnergyLight lamps producing light with all the intensity and quality of natural daylight. In fact, it can give you up to 10,000 lux, but without any unpleasant dazzle or glare. If you use the EnergyLight for just 30 minutes a day (there's no need to stare into the lamp all the time) at a light intensity of 10,000 lux, you'll notice an improvement in your energy level and mood after just 5 days' use.

Digital timer and display



The EnergyLight HF 3309 is equipped with a digital user interface with built-in user protocol for easy use, which also gives direct feedback on the amount of energy gained.

Light intensity dimmer



Fight winter blues and get the energy you have in summer while using the appliance at a convenient

level of light you need for easy-reading etc. The EnergyLight features a light intensity dimmer so it can easily be adjusted between zero to about 10,000 lux (at 15-20 cm) matching the intensity of your personal preference.

Adjustable angle



The appliance can be positioned at various angles for comfortable use.

UV-free



The EnergyLight uses the full natural daylight spectrum. It meets European safety standards and features a UV-filtering screen. The EnergyLight is therefore safe to use.

Supported by the Medical Device Directive (MDD)



This device has proven health benefits, as certified by the Medical Device Directive (MDD).

Specifications

Technical specifications

- Power: 72 W
- Cord length: 3 m
- Frequency: 50/60 Hz
- Insulation: Class II (double isolation)
- Voltage: 220/230
- Type of lamps: 2 PL-L 36W Philips EnergyLight
- Lifetime of lamps: 10,000
- Type of ballast: HF-R 232PLT (High Frequency)
- Ballast operating frequency: 42,000 kHz

Weight and dimensions

- Product dimensions: 33x57x18 (WxHxD) cm
- Product weight: 3.5 kg
- Qnt. on Euro pallet: 30 (UK 39) pcs
- F-box dimensions (= A-box): 36x62.2x21.5 (WxHxD) cm
- F-box weight (= A-box): 4.4 kg

Logistic data

- CTV code: 884330901000
- Country of origin: Germany

Safety

- UV-free: No UV comes through the filtering screen
- No flickering: due to high frequency ballast avoiding headaches
- Cenelec certified

Easy to use

- Digital user interface: for feedback on energy and time
- On product user protocol
- On/off button
- Energybar.: indicating how much energy is gained
- Digital timer
- Pause function

Feel energetic in 30 min a day

- Original Philips EnergyLight: lamps: 10,000 lux at 20 cm, 2,500 Lux at 60 cm

Comfortable light

- Light intensity dimmer

Easy positioning

- Adjustable angle: 20 degrees
- Compact size: 18x33 (DxW) cm

Satisfaction indicator

- 94% satisfied users
- Best buy: award by STIWA (2003) as successor of HF 3305

Medical appliance

- Proven health benefits



Issue date 2022-09-14

Version: 2.0.1

© 2022 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com