

PHILIPS

IPL hair removal
system

Lumea

For men

TT3000/00



Prevent hair regrowth

enjoy smooth skin everyday

Philips Lumea IPL hair removal system uses gentle pulses of light to stimulate the hair root into a resting phase. Hair sheds naturally and regrowth is inhibited. Using it every two weeks enables you to enjoy smooth skin everyday.

Effective hair regrowth control

- Body hair regrowth control

Intense Pulsed Light (IPL) technology

- Suitable for home use

Gentle application

- Hairs are eliminated before they appear

Application areas

- For chest, stomach, legs, underarms, back and pubic area

Suitable hair and skin types

- NOT suitable for red, grey, light-blond hair & skin type VI

Highlights

IPL technology



Philips Lumea works with an innovative light-based technology called "Intense Pulsed Light" (IPL) derived from technology used in professional beauty salons for hair removal. Philips Lumea now brings this innovative new technology, suitable for home use, to be used by you wherever and whenever you want.

Lamp does not need replacement



Philips Lumea comes equipped with a durable lamp with a lifetime of 50,000 flashes*

Suitable hair and skin types



Philips Lumea works effectively on (naturally) dark-blond, brown and black hair. Philips Lumea is NOT effective for (naturally) red, grey and light blond hair. Philips Lumea is NOT suitable for very dark skin (skin type VI)

Smoother results than a blade




Hairs are eliminated before they appear

Application body areas



Philips Lumea for Men has been developed in close cooperation with leading dermatologists for a safe, effective and gentle application on arms, underarms, shoulders, stomach, pubic area (excluding the scrotum), legs, chest and back

Body hair regrowth control

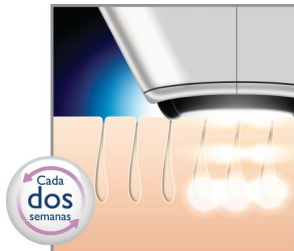

Philips Lumea for Men applies gentle pulses of light to the hair root. As a consequence the hair sheds naturally and hair regrowth is inhibited. Repeating this treatment every two weeks for the underarms, pubic area, chest, stomach, legs and back will keep your skin smooth every day. For a most effective treatment pre-shave the application areas with the included Bodygroom as long as you still experience regrowth

Specifications

Features

3cm² application window: For an effective treatment of large areas

5 light energy settings: Adjustable to your personal skin type

integrated UV-filter: Protects your skin from UV light

Ready-to-flash light: Indicates correct application on the skin

Accessories

Storage pouch: included

Cleaning cloth: included

AC/DC Adaptor/Charger: included

Technical specifications

Battery type: Lithium Ion

Frequency: 50-60 Hz

Power: 7.5 W

Voltage: 100-240 V

Light energy level: 2 - 6.5 J/cm² (depending on chosen setting)

Light spectrum: >570 nm

Logistic data

Country of origin: Hungary

