

PHILIPS

sonicare

Battery Sonicare
toothbrush

Essence

1 mode

1 brush head

HX5551/02



Better reach for a better clean

The only angled Sonic brush that gives you better reach for a better clean. Sonicare's unique dynamic cleaning action reaches deep between teeth and into hard-to-reach areas, while massaging gums for a brushing experience unlike any other.

Helps you meet recommended brushing guidelines

- Smartimer encourages 2 minute brushing time

Additional features

- Brush head holder
- Signals when it's time to recharge the battery

Better reach for a better clean

- Sophisticated brush head

Ease into the Sonicare experience

- Easy-start to acclimate to the Sonicare experience

Dynamic cleaning action

- Patented Sonic technology

Highlights


Patented Sonic technology

Sonicare's patented Sonic technology combines contour-fit bristles and wide brushing strokes with high-speed bristle motion (31,000 brush strokes per minute) to create a dynamic cleaning action. This gives you a brushing experience unlike any other. The Sonicare toothbrush is also clinically proven to help reduce gingivitis, improve gum health and significantly reduce coffee, tea and tobacco stains, for naturally whiter teeth.

Sophisticated brush head

Brush head with slim, angled neck and round-ended bristle tips, designed to fit the contours of the teeth for enhanced manoeuvrability.

Easy-start

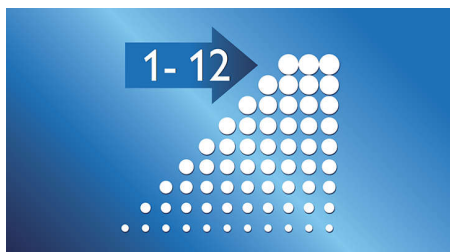

Gently increases power over first 12 uses, to ease into the Sonicare experience

Smartimer



Two-minute timer helps ensure Dental Professional recommended brushing time

Quadpacer



Quadpacer 30-second interval timer encourages even brushing

Recharge indicator



Signals when it's time to recharge the battery. Long-life battery provides up to two weeks of brushing, twice a day, between charges.

Brush head holder

Conveniently stores an extra brush head.



Specifications

Design and finishing

Colour: White

Power

Voltage: 110-220 V

Technical specifications

Operating time (full to empty): Up to 2 weeks

Battery: Rechargeable

Battery type: NiMH

Service

Warranty: 2-year limited warranty

Ease of use

Handle: Ergonomic design, Rubber grip for easy handling

Brush head system: Easy screw-on brush heads

Battery indicator: Light shows battery status

Brushing time: Up to 2 weeks

Items included

Brush heads: 1 e-Series standard

Charger: 1

Travel case: 1

Handles: 1 Essence

Cleaning performance

Health benefits: Helps improve gum health

Speed: Up to 62,000 brush movements/min,

Up to 31,000 brush movements/min

Performance: Up to 2 X more plaque removal*

Whitening benefits: Helps naturally remove stains

Timer: Quadpacer and Smartimer

Modes

Clean: For exceptional everyday cleaning

