



Philips Sonicare AirFloss
Electric Toothbrush &
Dental Flosser

w/5100 ProtectiveClean handle

1 nozzle
1 brush head



HX8272/21

Don't floss? Then AirFloss.

Removes plaque where brushing can't

For those who don't floss consistently, AirFloss is an easy way to start cleaning between teeth. AirFloss can be used with mouthwash or water and features unique air and micro-droplet technology to help remove plaque in hard-to-reach areas.

Improves oral health

- Helps prevent cavities between teeth
- Up to 99.9% plaque removal
- Cleans your entire mouth in just 30 seconds

Optimise your brushing

- Connects smart brush handle and smart brush heads
- Timers encourage thorough 2 minutes of brushing
- Whitens teeth in just one week
- Choose from three modes

Superior technology

- Air and micro-droplet technology


The easiest way to effectively clean between teeth

- An easy way to start a healthy habit

PHILIPS
sonicare

Highlights

Up to 99.9% Plaque Removal

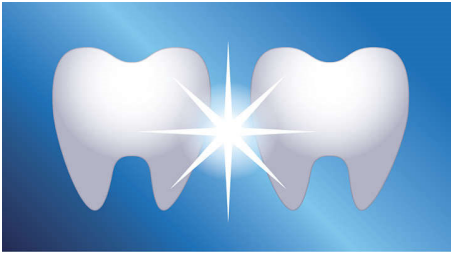

AirFloss uses microdroplets of air and water or mouthwash to remove up to 99% more plaque between teeth*

Helps prevent cavities



By gently bursting away plaque that brushing missed, Sonicare AirFloss helps prevent cavities from forming in the spaces between your teeth.

Micro-droplet technology



Our clinically proven results are possible thanks to our unique technology, which combines air and mouthwash or water to powerfully yet gently clean between teeth and along the gumline.

Easily start a healthy habit



Interproximal cleaning is very important for overall oral health. AirFloss is an easy way to clean deeper between teeth, helping to form a healthy habit.

30-Second Clean




With the AirFloss, cleaning your entire mouth takes just 30 seconds. Simply press the button once for a single burst, then move and repeat. To use the new auto-burst feature, hold down the button and move the guidance tip from space to space. It will automatically burst every second.

BrushSync technology



A microchip-enabled technology that detects and synchronises the smart brush head with the smart handle. The smart handle and smart brush head pair is a powerful combination that enables Smart replacement reminders.

Quadpacer and Smartimer

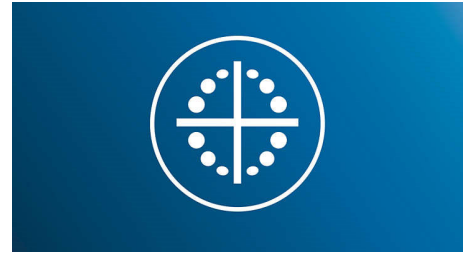


It takes just 2 minutes to thoroughly brush your teeth. Our Quadpacer lets you know when you've spent the optimal amount of time on each section of your mouth while our Smartimer signals when your total time is up. Together they help you reach your recommended brushing time, every time.

Naturally whiter teeth



Click on the W2 Optimal White brush head to remove surface stains and reveal a whiter smile. With its densely packed central stain-removal bristles, it's clinically proven to whiten teeth in just one week.

Three modes

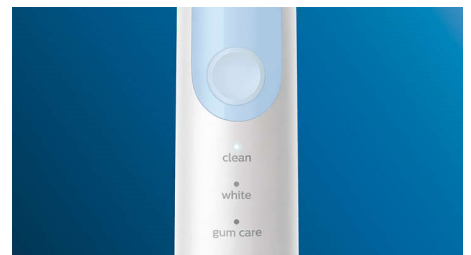

With a choice of three modes, this toothbrush allows you to customise brushing to your own needs. Clean mode is standard for superior cleaning. White is the ideal mode for removing surface stains. And Gum Care mode adds an extra minute of reduced-power brushing, so you can gently massage your gums.

Specifications

Design and finishing

- Colour: White and transparent green, White and Light Blue

Power

- Voltage: Multi-voltage charger

Technical specifications

- Operating time: Replace nozzle every 6 months
- Battery: Rechargeable
- Battery type: Lithium ION
- Charge time: 24 hours to full charge

Service

- Warranty: 2-year limited warranty

Ease of use

- Handle: Slim ergonomic design
- Battery indicator: Flashing indicates low battery

- Battery Life: Up to 2 weeks
- Nozzle attachment: Easily snaps on and off

Items included

- Charger: 1
- Brush heads: 1 W2 Optimal White
- AirFloss Handle: 1
- Handle: 1 ProtectiveClean
- AirFloss Nozzle: 1

Cleaning performance

- Auto-burst feature: Cleans entire mouth in 30 secs
- Speed: Up to 62,000 brush movements/min
- Whitening benefits: Whitens teeth in 1 week

Modes

- Clean: For exceptional everyday cleaning
- Gum Care: Gently massages gums
- White: Removes surface stains



Issue date 2023-02-15

Version: 5.0.2

EAN: 08 71010 38873 31

© 2023 Koninklijke Philips N.V.
All Rights reserved.

Specifications are subject to change without notice.
Trademarks are the property of Koninklijke Philips N.V.
or their respective owners.

www.philips.com

* Removes up to 7x more plaque vs. a manual toothbrush