

# PHILIPS

## sonicare

Bluetooth® connected  
toothbrush - Trial

FlexCare Platinum  
Connected

Built-in Bluetooth®

Smart brushing sensors

3 modes, 3 intensities

3 brush heads



HX9194/08



## Reveal their habits for improved brushing routines App turns brushing data into feedback and coaching

The new Philips Sonicare FlexCare Platinum Connected, powered by Smart Sensor Technology, reveals your patients' brushing habits and helps coach them into a better oral care routine.

### Provides a superior clean


- AdaptiveClean brush head
- BrushPacer and SmarTimer for thorough brushing
- Customizable modes and intensities
- Extended battery life and convenient travel case
- InterCare brush head

### Track your brushing with real-time guidance

- State-of-the-art 3D mouth-map
- Brush head monitor
- Goal-setting and progress reports
- Powered by Smart Sensor Technology
- Highlights trouble spots
- Custom guidance with Touch-Up

# Highlights

## Smart Sensor Technology

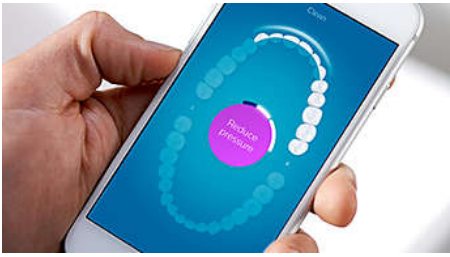

Three smart sensors send personalized brushing data to the Philips Sonicare app. Location sensor shows patients where they're brushing too little for better coverage. Scrubbing sensor helps guide patients to reduce scrubbing for optimized technique and a deeper, more gentle clean. Pressure sensor tracks brushing pressure and gently pulsates the handle to tell patients when to adopt a softer touch.

## State-of-the-art 3D mouth-map




A personalized 3D mouth-map provides post-brushing analysis so patients can track and refine their routine and identify trouble spots that require extra attention.

## Touch-Up




The Touch-Up feature displays customized guidance based on actual brushing data in order to immediately address missed areas.

## Highlights trouble spots

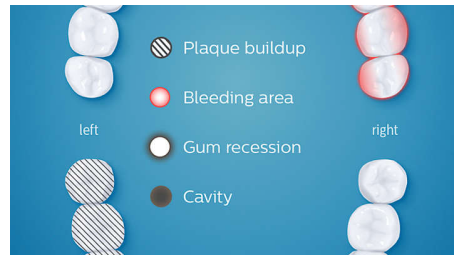

Based on advice from their dental professional, patients can highlight and address areas on the personalized mouth-map in need of extra care due to plaque buildup, gum recession or another issue.

## Brush head monitor

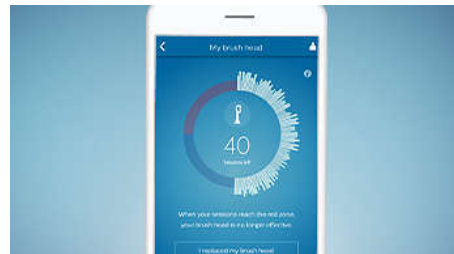

Using actual brushing time and pressure, the brush head monitor alerts patients in three months or less when it's time to replace their brush head, ensuring continued effectiveness for the best results possible.


## Goal-setting, progress reports



The Philips Sonicare app maintains a running history of a patient's personalized brushing data, allowing them to easily review performance and set daily and long-term goals. After 30 days, patients can generate an in-app progress report summarizing their

brushing performance, with the added option of sharing reports with their dental professional via email.

## BrushPacer and SmartTimer

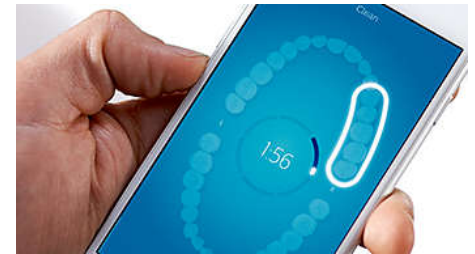

The BrushPacer tracks brushing through six segments of the mouth for a complete clean, while the SmartTimer ensures that patients brush for the full two-minutes as recommended by dental professionals.

## AdaptiveClean brush head



Soft, flexible sides made of rubber allow the bristles to conform to the unique shape of patients' teeth and gums to provide up to 4x more surface contact\*. The brush head's flexibility helps the bristles gently track along the gumline and absorb excessive brushing pressure. Enhanced sweeping motion delivers a unique in-mouth feel and superior cleaning along the gumline and in hard-to-reach areas\*. Proven to remove 10x more plaque than a manual toothbrush in hard-to-reach areas\*.

# Specifications

## Design and finishing

Color: White

## Connectivity

Bluetooth® wireless technology: Connected brushing app

## Power

Voltage: 110-220 V

## Technical specifications

Operating time (full to empty): Up to 2 weeks

Battery: Rechargeable

Battery type: Lithium ION

## Service

Warranty: 2-year limited warranty

## Compatibility

Android compatibility: Android phones, Bluetooth 4.0 enabled tablets

iOS compatibility: iPad 3rd Gen or higher, iPhone 4S or higher, with iOS7 operational system

## Ease of use

Battery indicator: Light shows battery status

Display: Illuminated display

Brush head system: Easy click-on brush heads

Handle: Slim ergonomic design

Brushing time: Up to 2 weeks

## Items included

Handles: 1 FlexCarePlatinum Connected

Charger: 1

Travel case: 1

Brush heads: 1 AdaptiveClean, 1 InterCare

standard, 1 ProResults gum health

Antibacterial tongue spray: 2oz

Tongue brush: 1

## Cleaning performance

Performance: Remove up to 10x more plaque\*\*

Health benefits: Improves gum health in only two weeks

Speed: Up to 62000 brush movement/min

Pressure sensor: Vibrates handle to alert user

Timer: BrushPacer and SmarTimer

Whitening benefits: Helps naturally remove stains

## Modes

Clean: For exceptional everyday clean

3 intensities: Low, Medium, High

White: Removes surface stains

Deep Clean: For an invigorating deep clean

