

PHILIPS

sonicare



Sonic electric
toothbrush

3 Series gum health

1 mode

3 intensity settings

2 brush heads



HX6632/18

Improves gum health up to 100%*

*more than manual brushing

Provides superior plaque removal and gum care in an easy-to-use power toothbrush to help prevent and reduce the early signs of gum disease while offering your most comfortable brushing experience

Guides you to follow dentist recommendations

- SmarTimer helps you meet recommended brushing time
- Quadpacer interval timer encourages thorough brushing

Improves gum health better vs. a manual toothbrush

- Better than a manual toothbrush for gingivitis
- Removes up to 6X more plaque along the gumline
- Sonicare dynamic cleaning action drives fluid between teeth

Offers a customized brushing experience

- Three intensity settings for a gentle experience

Recommended by dental professionals world wide

- Proven to improve gum health 3 ways
- Safe on orthodontics, dental work and veneers

3 Series gum health

Sonic electric toothbrush

HX6632/18

Highlights

Improves gum health 3 ways




Clinically proven to improve gum health three ways: reduces redness, swelling, bleeding (sites)

Prevents & reduces gingivitis



With superior plaque removal along the gum line, this toothbrush can help reduce gingivitis up to 100% better than a manual toothbrush in just two weeks

Three intensity settings

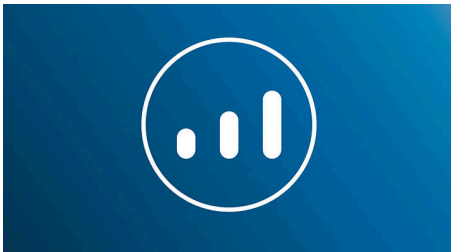

With three customizable intensity settings (low, medium & high) you can select the right power level to ensure that brushing along the gumline is comfortable. After 1.5 seconds of brushing, you will have to push twice on button to switch mode. Once to pause and second to switch.

Up to 6x better plaque removal



Removes up to 6X more plaque along the gumline than a manual toothbrush.

Sonic technology



Our sonic technology creates a unique dynamic cleaning action that drives fluid deep between teeth and along the gumline for a thorough-yet-gentle clean

Safe and gentle



The gentle motion of this Philips Sonicare power toothbrush is safe for dental work, including orthodontics, implants and veneers, so you can brush without worry.

SmarTimer



This Philips Sonicare power toothbrush features a special SmarTimer to help you meet the dentist-recommended brushing time of two minutes

Quadpacer



Features a Quadpacer interval timer to help encourage thorough brushing of the entire mouth. The toothbrush beeps every 30 seconds, signifying it's time to move to another area of the mouth.

3 Series gum health

Sonic electric toothbrush

HX6632/18

Specifications

Power

Voltage: 110-220 V

Technical specifications

Operating time (full to empty): Up to 3 weeks**

Battery: Rechargeable

Battery type: Lithium ION

Design and finishing

Color: Persian Green

Service

Warranty: 2-year limited warranty

Ease of use

Handle: Slim ergonomic design

Brush head system: Easy click-on brush heads

Battery indicator: Light shows battery status

Brushing time: Up to 3 weeks**

Items included

Travel case: 1

Brush heads: 2 ProResults gum health

Handles: 1 3 Series gum health

Charger: 1

Mouth rinse: 27 ml Listerine

Cleaning performance

Health benefits: Improves gum health

Performance: Removes up to 6X more plaque*

Speed: Up to 31000 brush movement/min

Timer: Quadpacer and SmarTimer

Modes

Clean: For exceptional everyday clean

3 intensities: Low, Medium, High

