

# PHILIPS



IPL - Hair removal  
device

Lumea Series 7000

For body

With skin tone sensor

15 minutes to treat lower legs



SC1994/80

## 12 months of hair-free smoother skin\*

Achieve 92% hair reduction in 3 treatments^

Get silky smooth skin from the comfort of your home with Philips Lumea Series 7000. Safe, effective, and easy to use - order yours today and enjoy long-lasting results.

### Long lasting hair removal with IPL technology

- Suitable for a wide variety of hair and skin types
- FREE Lumea App for a personal coach at your fingertips
- Developed with dermatologists to be safe, easy and effective.

### Safe, fast and effective on body

- Integrated skin tone sensor
- Achieve 92% hair reduction in as little as 3 treatments.

### Maximum freedom, minimum hassle

- Large treatment window for fast application


## Lumea Series 7000

IPL - Hair removal device

SC1994/80


# Highlights

### Developed with dermatologists

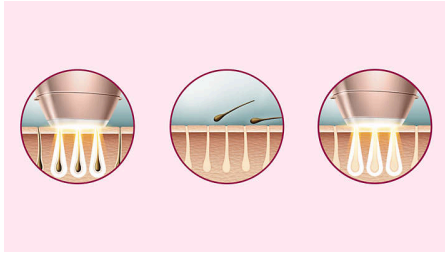

As a leader in health technology, Philips developed Lumea IPL in consultation with dermatologists, for easy and effective use in the safety of your home. Derived from technology used in professional salons, Lumea IPL offers gentle treatment, even on sensitive areas. The treatment leaves your skin beautifully hair free and touchably smooth.

### Proven effectiveness

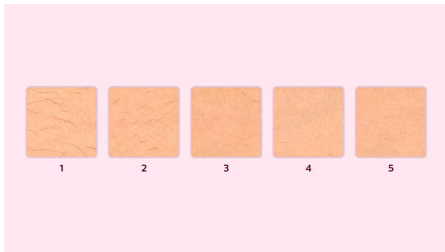

The first 4 treatments must be carried out every 2 weeks. Subsequently, to maintain your results, simply touch up every 4 weeks. With just 8 touch up treatments you can be hair-free for at least 3 months.

### Suitable hair and skin types



Philips Lumea Advanced works effectively, easily and safely on a wide variety of hair and skin types. It works on naturally dark blonde, brown and black colored hairs and on skin tones from very white to mid brown. As with other IPL based treatments, Philips Lumea cannot be used to treat white / grey, light blonde or red hair and is not suitable for very dark skin. This is due to the high contrast required between the pigment in hair color and pigment in the skin tone.


### Skin Tone Sensor

The integrated skin tone sensor measures the treated skin complexion at the beginning of each session and occasionally during the session. If it detects a skin tone that is too dark for this product it will automatically stop emitting pulses.

### Fast treatment on legs

Big treatment window for fast application on large areas like legs.

### FREE Lumea App

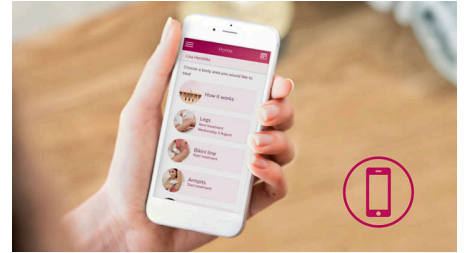

The Lumea App is free to download and offers a unique and personalized IPL experience – just for you. It is your very own personal coach to ensure you are getting the most out of your Lumea and using it in the correct way to achieve long-lasting results. The App helps you create your personalized treatment schedule per body area, with tips and advice during each treatment. All you need to do is check the app regularly to make sure you are up to date with any notifications or reminders.

# Lumea Series 7000

IPL - Hair removal device

SC1994/80

## Specifications

### Application areas

**Body areas:** Arms, Legs, Belly, Armpits, Bikini

### Safety and adjustable settings

**Integrated UV filter:** Protects skin from UV light

**5 light energy settings:** Adjustable to your skin type

**Integrated safety system:** Prevents unintentional flashing

**Skin Tone sensor:** Detects your skin tone

### Application mode

**Slide and Flash:** For easy motion application

**Stamp and Flash:** For treatment on curvy areas

**Corded / cordless use:** Corded

### Power

**Voltage:** 100-240 V

**Corded use**

### Consumer Trade Item

**Height:** 26.6 cm

**Width:** 19 cm

**Length:** 12 cm

**Net weight:** 0.826 kg

**Gross weight:** 1.08 kg

**GTIN:** 8710103982135

**Country of origin:** HU

**Harmonized system code:** 85437090

### Outer Carton

**Height:** 19.7 cm

**Width:** 24.8 cm

**Length:** 28 cm

**Net weight:** 1.652 kg

**Gross weight:** 2.402 kg

**GTIN:** 18710103982132

### Attachments

**Body attachment (4cm2):** For use below the neckline

### Service

**Warranty:** 2 years global warranty + 1 extra year of warranty upon product registration within 90 days

### Pallet Layer

**Height:** 19.7 cm

**Width:** 80 cm

**Length:** 120 cm

### Items included

**Instructions for use:** User manual

**Storage:** Pouch

**Adapter:** 24V / 1500mA

**Accessory:** Cleaning cloth

### Application time

**Bikini:** 4min

**Lower leg:** 15 min

**Armpit:** 2,5 min

