



Philips Brilliance  
LCD monitor with Clinical  
D-image

Clinical review display  
24.1" (61 cm)  
1920 x 1200 (2 MP)



C240P4QPYEW

## Clinical Review Display with Brilliant Images

This Philips display is ideal for everyday clinical care environments. Its D-image preset gives you consistent DICOM greyscale medical images. Its medical grade power supply ensures basic safety.

### Brilliant performance

- Built-in stereo speakers for multimedia
- DisplayPort offers audio and video over a single, long cable
- IPS technology for full colours and wide viewing angles
- USB Hub with four ports for easy connections
- PowerSensor saves up to 80% energy costs

### Clinical Review

- Design complements clinical environment
- Medical-grade power supply ensures safety with low emission
- Clinical D-image pre-set for DICOM part 14 Greyscale images
- SmartImage Clinic for common display settings

### Designed for people

- Low bezel-to-table height for maximum reading comfort
- SmartErgoBase enables user-friendly ergonomic adjustments

**PHILIPS**

LCD monitor with Clinical D-image  
Clinical review display 24.1" (61 cm), 1920 x 1200 (2 MP)

# Highlights

## Clinical D-image



Monitors must show medical images consistently with high quality to achieve reliable interpretations. The rendering of medical Greyscale images on standard monitors is mostly inconsistent at best, making them unsuitable for use in a clinical environment. Philips clinical review displays with clinical D-image pre-set are factory-calibrated to give DICOM part 14 compatible Greyscale standard display performance. By using high-quality LCD panels with LED technology, Philips offers you consistent and reliable performance at an affordable price. Further info on <http://medical.nema.org/>


## SmartImage Clinic

SmartImage is an exclusive, leading-edge Philips technology that analyses the content displayed on your screen, and based on different applications/images displayed on a monitor, SmartImage optimises display content by dynamically adjusting brightness, contrast, colour and sharpness in real time.

## Medical-grade power supply


Basic safety, essential performance and reliability are essential requirements in a professional clinical environment. Unlike standard displays, Philips clinical review displays employ special medical-grade power supplies in order to meet these requirements. It complies with the EN/IEC 60601-1 standard ensuring basic safety and essential performance. It also complies with the EN/IEC 60601-1-2 collateral standard for medical-grade electromagnetic tests and compatibility. Many countries, hospitals and clinics mandate compliance to these standards as a requirement.

## IPS technology

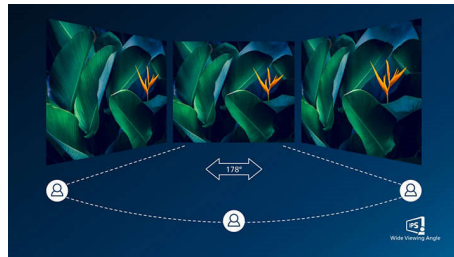

IPS displays use advanced technology that gives you extra-wide viewing angles of 178/178 degrees, making it possible to view the display from almost any angle — even in 90-degree Pivot mode! Unlike standard TN panels, IPS displays gives you remarkably crisp images with vivid colours, making it ideal not only for Photos, films and web browsing, but also for professional applications that demand colour accuracy and consistent brightness at all times.

## PowerSensor

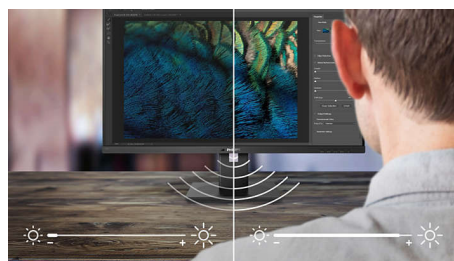

PowerSensor is a built-in 'people sensor' that transmits and receives harmless infrared signals to determine if the user is present and then automatically reduces monitor brightness when the user steps away from the desk, cutting energy costs by up to 80 percent and prolonging monitor life

## DisplayPort

DisplayPort is a digital link from PC to monitor without any conversion. With higher capabilities than DVI standard, it is fully capable to support up to 15 metre cables and 10.8 Gbps/sec data transfer. With this high performance and zero latency, you get the fastest imaging and refresh rates — making DisplayPort the best choice for not only general office or home use, but also for the more demanding gaming and films, video

C240P4QPYEW/00

editing and more. It also keeps interoperability in mind via the use of various adapters.

## SmartErgoBase



The SmartErgoBase is a monitor base that delivers ergonomic display comfort and provides cable management. The base can swivel, tilt and rotate to various angles to ensure maximum comfort. The height-adjustable stand guarantees the optimal viewing level, reducing the physical strains of a long workday, while cable management reduces cable clutter and keeps the workspace neat and professional.

## Four Ports USB 2.0 HUB

USB Hub allows user to conveniently connect their plug and play multimedia devices such as USB memory devices, camera, portable HDD, Web camera, PDA, Printer and many other devices which have USB connection available. The conveniently located USB 2.0 Hub on the monitor allows USB 2.0 signals to pass through to the computer. Note that many devices like cameras and HDD may need to be powered On independently, as they have higher power requirements than the monitor USB Hub can provide.

## Built-in stereo speakers



A pair of high-quality stereo speakers built into a display device. They may be visible front firing, or invisible down firing, top firing, rear firing etc. depending on model and design.



# Specifications

## Connectivity

- Signal Input: DisplayPort 1.2, DVI-D (digital, HDCP), VGA (Analogue)
- USB: USB 2.0 x 4
- Sync Input: Sync on Green, Separate Sync
- Audio (In/Out): PC audio-in, Headphone out
- Built-in speakers: 2 W x 2

## Picture/Display

- Panel Size: 24 inch/61 cm
- Aspect ratio: 16:10
- LCD panel type: IPS technology
- Backlight type: W-LED system
- Pixel pitch: 0.270 x 0.270 mm
- Optimum resolution: 1920 x 1200 @ 60 Hz
- Brightness: 300 cd/m<sup>2</sup>
- DICOM calibrated brightness: 265 cd/m<sup>2</sup>
- Display colours: 16.7 M
- Contrast ratio (typical): 1000:1
- Response time (typical): 14 ms
- Viewing angle: 178° (H)/178° (V), @ C/R >= 10
- Picture enhancement: SmartImage Clinic
- Effective viewing area: 518.4 (H) x 324 (V)
- Scanning Frequency: 24–94 kHz (H)/48–85 Hz (V)
- sRGB
- DICOM: Clinical D-image, DICOM part-14

## Convenience

- Plug and Play Compatibility: DDC/CI, Mac OS X, sRGB, Windows 8.1 / 8 / 7
- User convenience: SmartImage Clinic, PowerSensor, Menu, Power On/Off, 4:3/Wide
- OSD Languages: English, French, German, Italian, Portuguese, Russian, Simplified Chinese, Spanish, Brazil Portuguese, Czech, Dutch, Finnish, Greek, Hungarian, Japanese, Korean, Polish, Swedish, Traditional Chinese, Turkish, Ukrainian
- Other convenience: Kensington lock, VESA mount (100 x 100 mm)
- Control software: SmartControl

## Stand

- Height adjustment: 130 mm
- Pivot: 90 degree
- Swivel: -65/65 degree
- Tilt: -5/20 degree

## Power

- Off mode: Zero watts with Zero switch
- On mode: 31.4 W (typ.)
- Standby mode: 0.5 W (typ.)
- Power LED indicator: Operation - White, Standby mode - White (flashing)
- External power supply: 100-240VAC, 50-60 Hz (Medical cert.: IEC60601-1/EN60601-1)

## Dimensions

- Packaging in mm (W x H x D): 632 x 457 x 286 mm
- Product without stand (mm): 555 x 388 x 65 mm
- Product with stand (max height): 555 x 550 x 244 mm

## Weight

- Product with packaging (kg): 9.80 kg
- Product with stand (kg): 6.97 kg
- Product without stand (kg): 4.64 kg

## Operating conditions

- Altitude: Operation: +12,000 ft (3658 m), Non-operation: +40,000 ft (12,192 m)
- Temperature range (operation): 0°C to 40°C °C
- MTBF: 30,000 hour(s)
- Relative humidity: 20–80% %
- Temperature range (storage): -20°C to 60°C °C

## Sustainability

- Environmental and energy: PowerSensor, RoHS, WEEE
- Recyclable packaging material: 100 %
- Specific Substances: PVC/BFR free housing, Mercury free

## Compliance and standards

- EMC: IEC/EN60601-1-2, CE Mark
- Ergonomic/Others: JIS Z2801, TUV/GS, TUV Ergo, ISO 13485 certified, CCC, CECP, RCM, UL/cUL, TCO certified
- Safety: IEC/EN60601-1, CE Mark

## Cabinet

- Finish: Texture
- Foot: White
- Front bezel: White
- Rear cover: White



Issue date 2024-04-24

Version: 5.0.1

EAN: 87 12581 68599 7

© 2024 Koninklijke Philips N.V.  
All Rights reserved.

Specifications are subject to change without notice.  
Trademarks are the property of Koninklijke Philips N.V.  
or their respective owners.

[www.philips.com](http://www.philips.com)

\* Please note: This display is not certified to be used in patient vicinity.  
Any use for primary diagnostics is the sole responsibility of the hospital or clinic operator only.  
\* DICOM is the registered trademark of the NEMA for its standards publications relating to digital communications of medical information.  
\* "IPS" word mark / trademark and related patents on technologies belong to their respective owners.