

# PHILIPS

## sonicare



Oral Irrigator

Cordless Power  
Flosser 3000

Quad Stream technology

Pulse waves guide as you go



HX3826/33

## Flossing reinvented with QuadStream technology

Removes up to 99.9% of plaque\* in treated areas

Effortlessly thorough cleaning between teeth to improve your gum health. Unique Quad Stream tip cleans more area with less effort for more effective flossing every time, so it's easy to get flossing right.

### Easy and convenient

- Fast universal charging with USB-A to Small Plug cable
- Easy-to-fill 250 ml reservoir
- Effortless technique
- Standard tip for spot cleaning

### Faster, more effective flossing

- 2 flossing modes, 3 intensities
- Complete clean in 60 seconds
- Pulse Wave technology guides from tooth to tooth
- Quad Stream technology for faster, more effective flossing

### An effortless clean with advanced technology

- Quadstream technology: removes plaque from below the gumline

# Cordless Power Flosser 3000

Oral Irrigator

HX3826/33

## Highlights

### Quad Stream technology




Unique X-shaped Quad Stream tip separates flow into 4 water jets, which cover more area between teeth and along the gum line. Get a faster, deeper clean without the hassle of string floss.

### Cleans from below the gumline



The Quadstream technology removes up to 99% of plaque from 6mm deep periodontal pockets.

### Pulse Wave technology

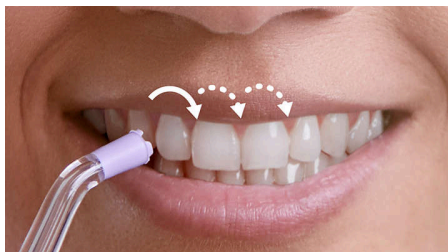

Gentle pulses of water keep you going and guide you from tooth to tooth in Deep Clean mode — so you get it right every time.

### 2 modes, 3 intensities



Choose the clean that suits your smile. Clean mode delivers continuous flow for effective everyday cleaning. Deep Clean pulses for more thorough cleaning. Adjustable intensity level adapts for comfort.

### Complete clean in 60 seconds




Access those hard-to-reach areas with a nozzle that rotates 360 degrees, so that you can use the device in any orientation. Get a thorough 360 clean in just 60 seconds from start to finish.

### Fast charging with USB-A to Small Plug cable




Convenient charging with USB-A to Small Plug cable. You can clean for up to 14 days on just one charge.

### 250ml reservoir

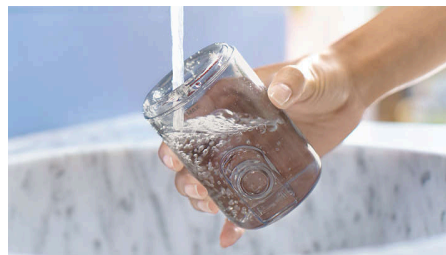

Reservoir capacity of 250 ml holds enough water for a recommended 60-second clean. Simply twist and remove or use the side-fill door for easy refilling.

### Effortless technique

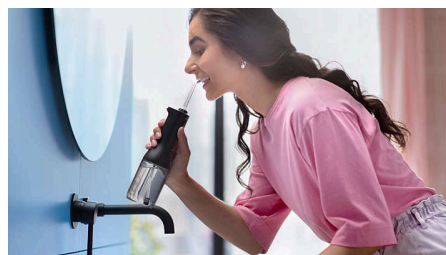

Just point, press and clean for more effective flossing. Aim the tip between your teeth and along the gumline and let it do all the work for you!

### Standard tip for spot cleaning



A standard single-stream tip boosts pressure for spot cleaning and removing food debris when you need it.

# Cordless Power Flosser 3000

Oral Irrigator

HX3826/33

## Specifications

### Power

**Voltage:** Multi-voltage charger

### Technical specifications

**Energy Consumption:** Off mode <0.3 W, Standby without display <0.5 W

### Consumer Trade Item

**Height:** 24.7 cm

**Width:** 13.4 cm

**Length:** 9.2 cm

**Net weight:** 0.458 kg

**Gross weight:** 0.676 kg

**GTIN:** 8720689011150

**Country of origin:** CN

**Harmonized system code:** 96032100

### Outer Carton

**Height:** 27.2 cm

**Width:** 15.1 cm

**Length:** 29.9 cm

**Net weight:** 1.374 kg

**Gross weight:** 2.328 kg

**GTIN:** 18720689011157

### Pallet

**Height:** 96.6 cm

**Width:** 80 cm

**Length:** 120 cm

**Net weight:** 82.44 kg

**Gross weight:** 164.94 kg

**GTIN:** 68720689011152

### Service

**Warranty:** 2-year limited warranty

### Pallet Layer

**Height:** 27.2 cm

**Width:** 80 cm

**Length:** 120 cm

### Ease of use

**Handle:** Slim ergonomic design

**Nozzle attachment:** Easily clicks on and off

### Performance

**Cleaning:** Clean entire mouth in 60-90sec

### Items included

**Cordless Power Flosser:** 1

**F1 Standard nozzle:** 1

**F3 Quad Stream nozzle:** 1

**Travel bag:** 1

**USB charging cable:** 1

**USB wall adaptor:** 1

**USB wall adaptors:** 1

### Modes

**Clean:** For a thorough everyday clean

**Deep Clean:** For a deep interdental clean

**Intensities:** 3

**Deep Clean:** For a deep interdental clean

### Health benefits

**Plaque removal:** Removes up to 99% of plaque in treated areas\*

**Sub-gingival:** Removes up to 99% of plaque from below the gumline\*\*

### Design and finishing

**Color:** Black

