

PHILIPS

sonicare



Oral Irrigator

Cordless Power
Flosser 3000

Quad Stream technology

Pulse waves guide as you go



HX3826/31

Flossing reinvented with QuadStream technology

Removes up to 99.9% of plaque* in treated areas

Effortlessly thorough cleaning between teeth to improve your gum health. Unique Quad Stream tip cleans more area with less effort for more effective flossing every time, so it's easy to get flossing right.

Easy and convenient

- Fast universal charging with USB-A to Small Plug cable
- Easy-to-fill 250 ml reservoir
- Effortless technique
- Standard tip for spot cleaning

Faster, more effective flossing

- 2 flossing modes, 3 intensities
- Complete clean in 60 seconds
- Pulse Wave technology guides from tooth to tooth
- Quad Stream technology for faster, more effective flossing

An effortless clean with advanced technology

- Quadstream technology: removes plaque from below the gumline

Cordless Power Flosser 3000

Oral Irrigator

HX3826/31

Highlights

Quad Stream technology




Unique X-shaped Quad Stream tip separates flow into 4 water jets, which cover more area between teeth and along the gum line. Get a faster, deeper clean without the hassle of string floss.

Cleans from below the gumline



The Quadstream technology removes up to 99% of plaque from 6mm deep periodontal pockets.

Pulse Wave technology

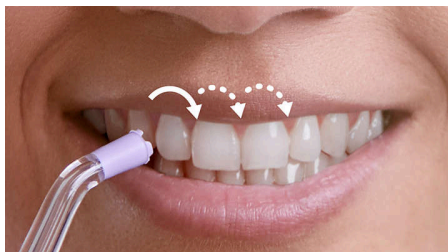


Gentle pulses of water keep you going and guide you from tooth to tooth in Deep Clean mode — so you get it right every time.

2 modes, 3 intensities



Choose the clean that suits your smile. Clean mode delivers continuous flow for effective everyday cleaning. Deep Clean pulses for more thorough cleaning. Adjustable intensity level adapts for comfort.

Complete clean in 60 seconds

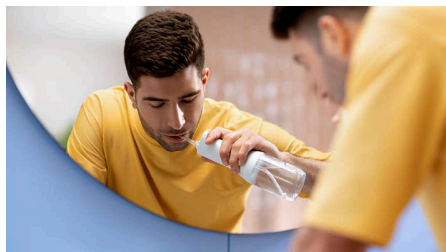


Access those hard-to-reach areas with a nozzle that rotates 360 degrees, so that you can use the device in any orientation. Get a thorough 360° clean in just 60 seconds from start to finish.

Fast charging with USB-A to Small Plug cable




Convenient charging with USB-A to Small Plug cable. You can clean for up to 14 days on just one charge.

250ml reservoir

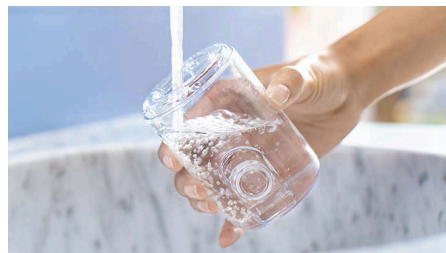


Reservoir capacity of 250 ml holds enough water for a recommended 60-second clean. Simply twist and remove or use the side-fill door for easy refilling.

Effortless technique

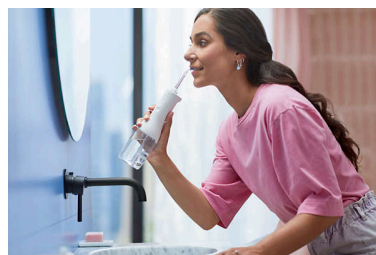

Just point, press and clean for more effective flossing. Aim the tip between your teeth and along the gumline — and let it do all the work for you!

Standard tip for spot cleaning

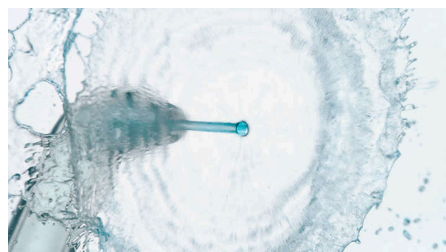

A standard single-stream tip boosts pressure for spot cleaning and removing food debris when you need it.

Cordless Power Flosser 3000

Oral Irrigator

HX3826/31

Specifications

Power

Voltage: Multi-voltage charger

Technical specifications

Energy Consumption: Off mode <0.3 W, Standby without display <0.5 W

Consumer Trade Item

Height: 24.7 cm

Width: 13.4 cm

Length: 9.2 cm

Net weight: 0.458 kg

Gross weight: 0.676 kg

GTIN: 8720689011129

Country of origin: CN

Harmonized system code: 96032100

Outer Carton

Height: 27.2 cm

Width: 15.1 cm

Length: 29.9 cm

Net weight: 1.374 kg

Gross weight: 2.328 kg

GTIN: 18720689011126

Pallet

Height: 96.6 cm

Width: 80 cm

Length: 120 cm

Net weight: 82.44 kg

Gross weight: 164.94 kg

GTIN: 68720689011121

Service

Warranty: 2-year limited warranty

Pallet Layer

Height: 27.2 cm

Width: 80 cm

Length: 120 cm

Ease of use

Handle: Slim ergonomic design

Nozzle attachment: Easily clicks on and off

Performance

Cleaning: Clean entire mouth in 60-90sec

Items included

F1 Standard nozzle: 1

F3 Quad Stream nozzle: 1

Travel bag: 1

USB charging cable: 1

USB wall adaptor: 1

Sonicare Cordless Power Flosser: 1

Modes

Clean: For a thorough everyday clean

Deep Clean: For a deep interdental clean

Intensities: 3 levels

Deep Clean: For a deep interdental clean

Health benefits

Sub-gingival: Removes up to 99% of plaque from below the gumline**

Design and finishing

Color: White

