

# PHILIPS

## sonicare



### Rechargeable Sonic Toothbrush

6500

Up to 10 x more plaque removal

3 modes, 3 intensities

Visible pressure alert

Philips Sonicare app



HX7411/01

## Gentle clean, even for sensitive teeth and gums

Up to 10 x more plaque removal\*

The Philips Sonicare 6500 features Next-Generation Sonicare technology, a dedicated Sensitive mode, 3 intensity levels, pressure alert and two brush heads: S2 Sensitive and W2 Whitening

#### Brushing made effortless

- Personalise your clean with 9 brushing settings

#### Brushing made effortless

- Elevate your routine with the Philips Sonicare app
- Protect your gums with the visible pressure alert
- Always know when it's time for a new brush head

#### Superior clean

- Next-Generation Sonicare Technology

#### Gentle, superior cleaning

- Gentle cleaning and whitening for sensitive teeth and gums

#### Designed for your lifestyle

- Make travel easy with our compact, ultra-durable travel case
- Optimise your brushing with the SmarTimer and BrushPacer
- Laser-welded for strength and durability


6500

Rechargeable Sonic Toothbrush

HX7411/01


# Highlights

## S2 Sensitive & W2 Whitening

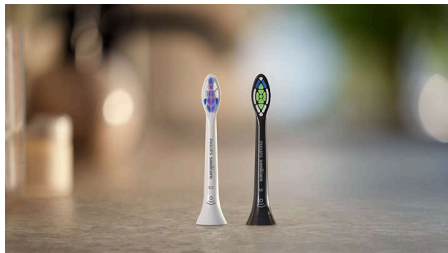

The Philips Sonicare Series 6500 electric toothbrush delivers a gentle brushing experience, even for sensitive teeth and gums. Thanks to the new S2 Sensitive Brush Head made using Sonicare's softest bristles, it is uniquely designed to provide a gentle clean that removes up to 10 x more plaque in hard-to-reach areas.\*\* In addition, the W2 Optimal White Brush Head makes it easy to maintain brightness between professional whitening treatments, while removing up to 7 x more plaque vs. a manual toothbrush.

## Next-Gen Sonicare technology

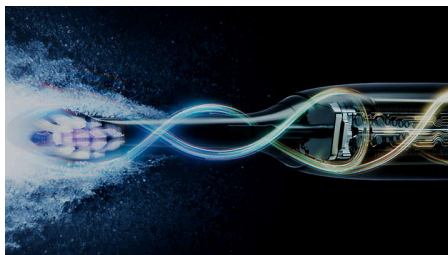


Enjoy a delightful brushing experience and new level of oral care with Next-Generation Sonicare technology. While you brush, the electric toothbrush's innovative magnetic system adjusts brushing power during each cleaning session, delivering optimal performance of precisely 62,000 bristle movements per minute, even in hard-to-reach areas. Meanwhile, dynamic fluid action supports the brush head's cleaning performance by gently rinsing along the gumline and between every tooth, helping you achieve a complete clean every time you brush.

## 9 brushing settings




Designed for sensitive teeth and gums, the Philips Sonicare Series 6500 electric toothbrush allows you to choose between 3 brushing modes, ensuring a clean that's right for you. Sensitive mode provides gentle, thorough cleaning, while Clean mode focuses on plaque removal for a superior clean. White mode works a bit harder to remove surface stains and polish your teeth for a whiter smile. For a personalised clean, you can combine any mode with one of 3 intensity levels, from High to Low, for a brushing session that's right for you.

## Light-ring pressure alert

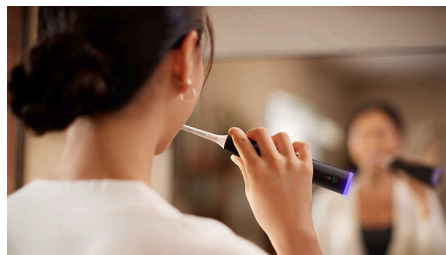

You may not notice when you're brushing too hard, but this electric toothbrush will. The smart optic pressure sensor precisely detects excessive brushing pressure, then flashes purple while reducing vibration. You'll always know when to adjust your brushing, and your gums will stay protected.

## Philips Sonicare app

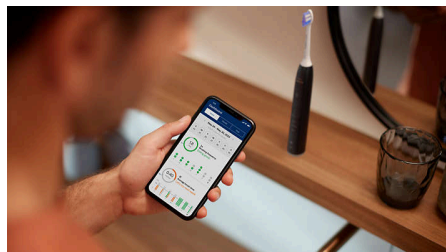

The Philips Sonicare app seamlessly pairs with your electric toothbrush and supports you in improving and maintaining your oral health. Our

app keeps track of your brushing progress, providing personalised guidance and progress reports so you can gain valuable insights about your brushing.

## 2-minute timer and brushing pacer



Achieve a complete clean in just two minutes with our SmarTimer, which ensures you brush for the dentist-recommended time. For an even clean throughout your mouth, the BrushPacer divides your session into segments, with a slight pause in vibration as your cue to move on. This supports you in achieving an even clean every time you brush.

## BrushSync replacement reminder



Replacing your electric toothbrush head regularly is an important part of maintaining good oral health. With BrushSync technology, it's easy to keep your brush head performing at its best. By tracking your usage frequency and brushing intensity, the innovative technology then reminds you when it's time for a new brush head.


6500

Rechargeable Sonic Toothbrush

HX7411/01

# Highlights

## Designed to last

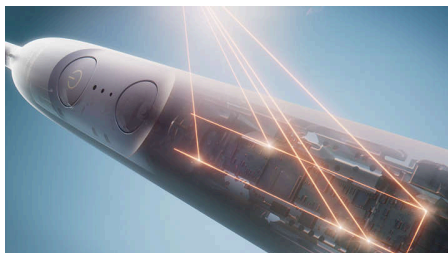

We've re-designed our Philips Sonicare electric toothbrushes from the inside out to be strong, lightweight and long-lasting. Thanks to our advanced laser-welding process, our electric toothbrushes are designed to brush with you for years to come.

## Up to 3-week battery life



This electric toothbrush features a high-performance long-lasting rechargeable battery designed to enhance your brushing experience. With a single charge, you can enjoy up to three weeks of regular brushing before needing to recharge, ensuring convenient and hassle-free oral care

## Compact travel case



Our redesigned travel case is compact enough to slip into a bag and strong enough to keep your Philips Sonicare electric toothbrush protected along the way. Not only does it keep your toothbrush protected, but it also conveniently holds two brush heads. Now you can travel confidently wherever your journey takes you.

6500

Rechargeable Sonic Toothbrush

HX7411/01

# Specifications

## Intensities

**High:** To boost your clean

**Medium:** For an everyday clean

**Low:** For sensitive teeth and gums

## Product Update

**Support Period:** Philips will provide the necessary security updates until 31 December 2026.

## Power

**Voltage:** 100-240 V

## Technical specifications

**Operating time (full to empty):** 21 days

**Battery:** Rechargeable

**Battery type:** Lithium ION

**Energy consumption:** Standby without display  
<0.5 W

## Design and finishing

**Colour:** Black

## Service

**Warranty:** 2-year limited warranty

## Compatibility

**Bluetooth wireless technology:** Connected brushing app

**Android compatibility:** Android phones with Bluetooth 4.0 and above

**iOS compatibility:** iPhone with Bluetooth 4.0 and above

## Ease of use

**Handle compatibility:** Easy click-on brush heads

**Handle:** Slim ergonomic design

**Battery indicator:** Illuminated icon indicates battery life

**Timer:** SmarTimer and BrushPacer for complete and thorough brushing

**USB-A Charging:** Power adapter not included

## Items included

**Travel case:** 1 slim travel case

**Charger:** 1

**Handle:** 1 6500 rechargeable toothbrush

**Brush head:** 1 S2 Sensitive, 1 W2 Optimal White

## Cleaning performance

**Plaque removal:** Up to 10 x more effective\*

**Speed:** 62,000 brush movements/min

**Gum Health:** Up to 5 x healthier gums\*\*

**Whitening:** Up to 2 x whiter teeth\*\*\*

## Modes

**Clean:** For exceptional everyday cleaning

**White:** To remove surface stains

**Sensitive:** For sensitive teeth and gums

## Smart sensor technology

**Pressure feedback:** Light ring illuminates purple, Reduced vibration and pulsing sound

**Replacement Reminder:** BrushSync lets you know when to replace your brush head

