

# PHILIPS

## AVENT



### Double Electric Starter Set

#### Hands-free

Mimics baby's drinking rhythm

Freedom while expressing

Easy to use, clean and assemble

Inclusive, seamless fit



SCD553/11

## Free your hands. Free your time.

A perfect fit for 99% of mums

Prepare for the new arrival with this Hands-free Double Electric Breast Pump Starter Set. It has everything you need to pump, store breast milk and feed baby from day one including breast shields and inserts, pads, bottles and storage bags.

#### Other benefits

- Everything you need to store, feed and more

#### Physical freedom while expressing

- Move freely with the ultra-light collection cups
- Wearable comfort for maximum movement

#### Mimics baby's drinking rhythm

- Express effectively\* with a rhythm that mimics your baby's
- Quiet hospital-strength pumping

#### Easy to use, simple to clean and assemble

- Track milk flow with the fully transparent collection cups
- Personalise your pumping session
- Easy to clean and assemble with few parts

#### Inclusive, seamless fit

- The complete range of breast shields and insert sizes

## Hands-free

Double Electric Starter Set

SCD553/11

# Highlights

### Mimics baby's drinking rhythm



Babies know how to drink best! That's why our hands-free electric breast pump mimics their natural drinking rhythm – so you can find your own rhythm for effective\* expression and optimal milk output. Pumping up to 85 times per minute, the Philips Avent Hands-free Breast Pump is 2 x faster than most other pumps\*\*.

### Hospital-strength motor




The nearly silent motor unit delivers the hospital-strength power you expect from a traditional breast pump, yet from a compact and wearable unit. With its rechargeable battery, the cordless motor unit pumps powerfully and discreetly anywhere.

### Unique SkinSense shields




Our unique SkinSense silicone breast shields are designed for ultimate comfort, softly shaping around your breasts using your body's natural warmth. The breast shields are made with food-grade silicone for complete peace of mind.

### Fits 99% of all mums

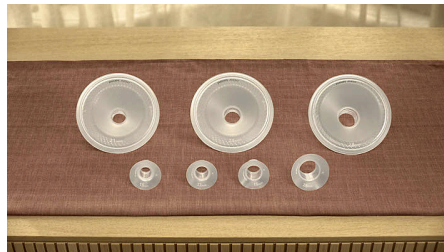

We support every mum to feel their best and having the right fit is an essential part of that. This starter set includes our full range of breast shields and shield insert sizes, fitting 99% of mums for ultimate comfort while pumping.

### Transparent collection

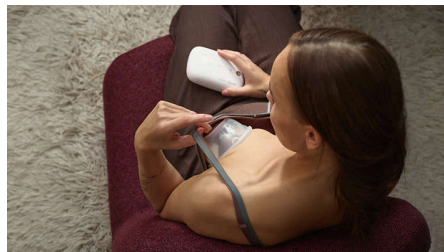

The see-through collection cups provide an unobstructed view inside. They help you correctly place your nipples for your best milk

output and visually track the volume of breast milk. For confident pumping from the first drop until finish.

### Ultra-light collection cups



Our collection cup weighs just a little over 100 grams when empty - up to 3 x lighter than a fully wearable breast pump. The cup sits comfortably in-bra with barely a sign that it's there.

### Wear it your way



Our pump comes with a braided carrying cord and apparel clip for wearing the motor unit cross-body or attached to your waistband. With the freedom to wear it your way, you stay in your flow knowing control is right by your side. Use the tube coupler and cord adjuster to precisely fine-tune your fit.


## Hands-free

Double Electric Starter Set

SCD553/11

# Highlights

### Personalised pumping

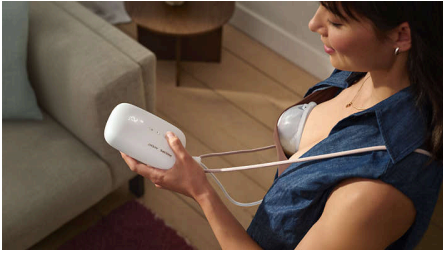

Get the most from every session with 8 suction levels for stimulating the flow of milk and 16 suction levels for expressing it. You can fine-tune as you pump and with your ideal setting automatically memorised, your pump is always ready to go.

### Easy to clean and assemble



We designed our pump and accessories to be as simple as possible. With only a few connecting pieces, each collection cup is a breeze to take apart, clean (by hand or in a dishwasher), disinfect and put back together.

### Starter set essentials



Experience convenience and flexibility with the starter set's feeding and storage products, designed to fit seamlessly into your routines. For extra comfort, the set also includes five pairs of breast pads.

## Specifications

### Power

**Source:** Rechargeable battery / USB-C (adapter not included)

### Material

**Cup, front cap, breast shield cover:**

Polypropylene

**Breast shield, membrane, valve, insert:** Liquid silicone

**Silicone tube:** Silicone

### What is included

**Disposable breast pad:** 10 pcs

**Breast milk bags:** 25 pcs

**USB-C cable:** 1 pcs

**Pouch:** 1 pcs

**Carrying cord:** 1 pcs

**Tube coupler:** 1 pcs

**Rechargeable motor unit (with a detachable clip):** 1 pcs

**Collection cup (assembled):** 2 pcs

**Silicone tubes:** 2 pcs

**Breast shield (dust) cover:** 2 pcs

**125 ml (4 oz) Natural Response bottle (with Flow 2 teat):** 1 pcs

**Nipple sizing card:** 1 pcs

**Inserts (15 mm, 17 mm, 19 mm, 26 mm):** 2 pcs each

**Breast shield (21 mm, 24 mm, 28 mm):** 2 pcs each

### Ease of use

**Breast pump use:** Easy to use, simple to clean and assemble

### Functions

**Flexibility in movement:** Hands-free pumping

**Personalisation settings:** 16 suction levels in expression mode, 8 suction levels in stimulation mode

**Motor unit:** Rechargeable battery / USB-C, Fully charging the battery takes up to 1.5 hours and lasts approx. 5 pumping sessions

**Breast shield and insert specifics:** Several sizes available for the best fit

