


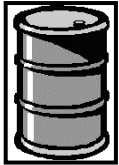


Product name:	HeartStart MRx Defibrillator	
Identification code	M3535A, M3536A, or M3536J	
Total weight (in Kg)	Up to 8 kgs	
Producer /Manufacturer	Name company:	Philips Medical Systems
	Address:	3000 Minuteman Road, Andover, MA
	Zip code:	01810
	Country:	USA
	Electronic info:	www.philips.pms.recycleinfo.com

Recycle info	Substances:	Location
Special attention	Item	Location
	<p>WARNING: Remove all power sources (AC, external batteries, or DC). See Figure 1. Wait at least 2 hours for any charge on the high voltage defibrillation capacitor and/or power supply capacitors to discharge before beginning disassembly. After waiting 2 hours, fully discharge the high voltage capacitor prior to touching inside the device. Dangerous voltages may be present on components and connections. Using an insulated probe, short between the red and white wire attachment points on the pc board. See Figure 2. Failure to do so could result in serious injury or death.</p>	Figure 1 & 2
Batteries  To be removed	Type:	Location
	External Lithium Ion (up to 2 batteries may be present)	Figure 1
	Internal Lithium	Figure 3
Hazardous  To be removed	LCD display backlight tubes contain mercury (Hg)	Figure 4 & 5
	Magnesium Casting	Figure 5
Fluids / Gases 	Item	Location
	None	

Locations as mentioned in the Passport


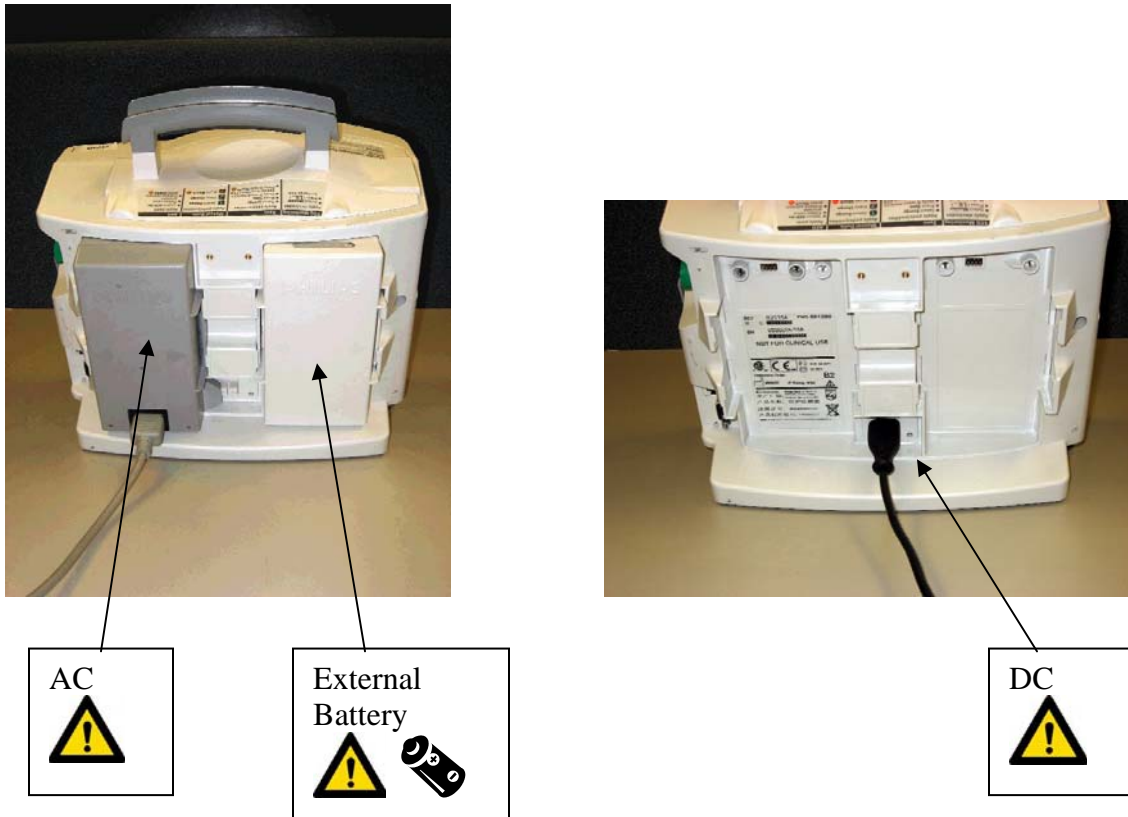


Figure 1

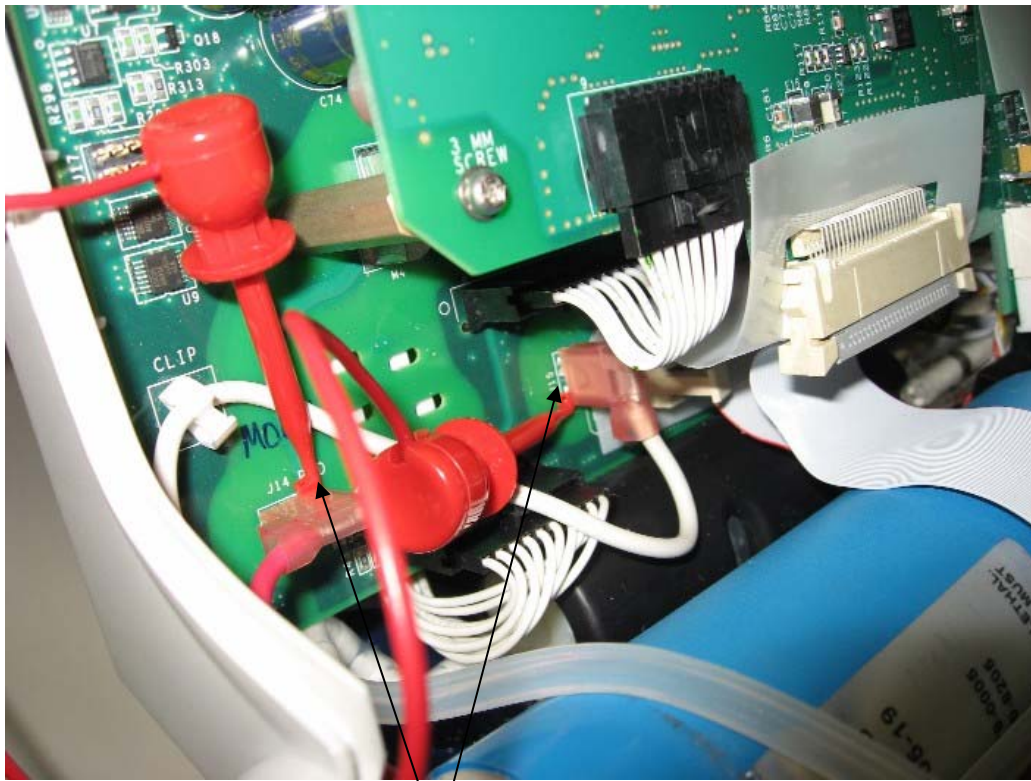

Discharge High
Voltage Capacitor



Figure 2



Internal
Lithium
Battery



Figure 3



LCD Display
(contains Hg)



Figure 4



Magnesium casting
and LCD Display are
located under pc board.



Magnesium casting



LCD Display
(contains Hg)



Figure 5