


<b>Product name</b>	<b>ST SENSE Breast Dx Coil</b>
Identification code(s)	9896-032-06151 ST SENSE Breast Dx 9896-032-06431 Phantom Kit, ST SENSE Breast Dx 9896-032-06421 Ancillary Kit, ST SENSE Breast Dx 4535-302-68731 FRU, ST SENSE Breast Dx 4535-302-68781 FRU, Phantom Bottle, ST SENSE Breast Dx 4535-302-68761 FRU, Cable Assembly, ST SENSE Breast Dx
Total NET weights [kg]	-06151: 5.7Kg, -06431: 11Kg, -06421:4.4Kg, -68731: 5.7Kg, -68781:5.7Kg, -68761: .66Kg
<b>Manufacturer</b>	<b>Invivo, Corporation. USA</b>
Address	3545 SW 47th Avenue
Zip code	Gainesville, FL 32608
Country	USA
Electronic info	Recycling passport ST SENSE Breast Dx Coil.pdf


### Recycling information

Special Attention	Items	Location
	None	n/a

Batteries	Items, type	Location
	None	n/a

Hazardous substances	Items	Location
	Several printed circuit boards inside the coil	Figure 1
	Beryllium copper contact pins on connectors	Figure 2
	Printed circuit board on cable assembly	Figure 3

Name: <b>Recycling passport ST SENSE Breast Dx Coil</b>	12NC: <b>9896 032 06151</b>	Version: <b>02</b> Doc.Date: <b>2008-11-06</b>
Author: <b>Faiz Ikramulla</b>	Page <b>1</b> of <b>3</b>	Group No: <b>5790</b> Section: <b>-</b> Status: <b>Authorized</b>

Fluids and gases	Items	Location
	Phantom Kit: , ST SENSE Breast Dx: 3 bottles NiCl	Figure 4
	Phantom Bottle:  NiCl: Water >1% Nickel Chloride hexahydrate <0.5% Sodium Chloride <0.5% (Weight %)	Figure 5

**Locations as mentioned in the Passport (Information in pictures)**

*Place sketch (es) or photo(s) here, and indicate issues using icons:*



Figure 1



Figure 2



Figure 3

Name: <b>Recycling passport ST SENSE Breast Dx Coil</b>	12NC: <b>9896 032 06151</b>	Version: <b>02</b> Doc.Date: <b>2008-11-06</b>
Author: <b>Faiz Ikramulla</b>	Page <b>2</b> of <b>3</b>	Group No: <b>5790</b> Section: <b>-</b> Status: <b>Authorized</b>

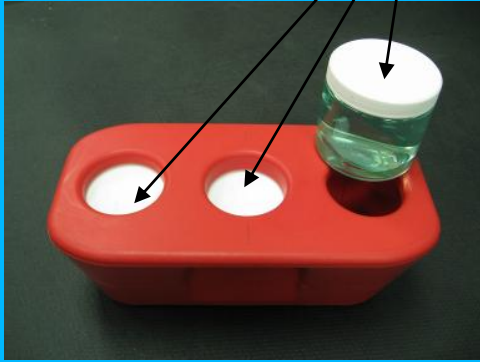


Figure 4

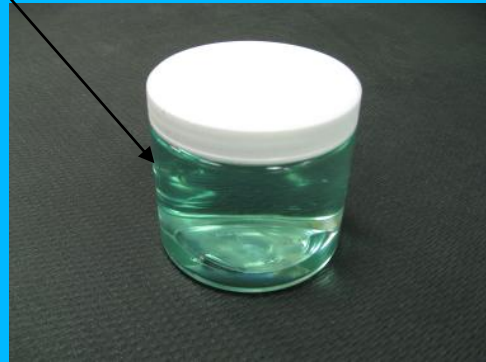


Figure 5

End of document.

Name: <b>Recycling passport ST SENSE Breast Dx Coil</b>	12NC: <b>9896 032 06151</b>	Version: <b>02</b> Doc.Date: <b>2008-11-06</b>
Author: <b>Faiz Ikramulla</b>	Page <b>3</b> of <b>3</b>	Group No: <b>5790</b> Section: -
		Status: <b>Authorized</b>