



Built for peace of mind

IntelliVue MX550 patient monitor

PHILIPS



Addressing the changing needs of healthcare

Philips IntelliVue patient monitors are designed to deliver accurate, real-time physiological patient data clinicians the world over have come to rely on when making decisions. Healthcare is a dynamic, demanding environment with constantly evolving needs. The MX series from Philips best-selling IntelliVue family aims to answer these challenges. You can now take advantage of efficient interfacing of bedside device data to electronic medical records (EMRs) – even during transport. And with an intuitive view of all available patient information, you can diagnose with confidence.

The consistent IntelliVue interface makes training and use as simple as possible. What's more, the IntelliVue family of patient monitors uses standard protocols to interface to other vendors' applications and work easily within your institution's existing IT infrastructure.



One intuitive view at the point of care


The IntelliVue MX550 pairs powerful bedside monitoring with the reassurance that a battery backup can bring. By supplying comprehensive patient information at a glance, the IntelliVue MX550 can make a real difference when multiple patients and priorities demand attention.

Monitoring needs change with patient acuity. That's why our devices are designed to monitor a wide range of vital signs in the ICU, on a sub-acute floor, or at the bedside. The IntelliVue MX550 establishes device interfacing and connectivity to the EMR at virtually any point of care in the hospital – even during transport or power loss.

Along with these exciting features, the IntelliVue MX550 offers the same high standard of patient monitoring for which our IntelliVue family is known.

See it clearly and quickly

Find what you need right away on the large 15"-wide touchscreen. You'll recognize the familiar, easy-to-use interface from your existing IntelliVue monitors, so you can spend more time providing care and less time on device training.

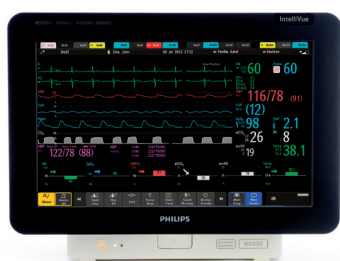

Whether moving a patient throughout the hospital or stationary at the bedside, you won't have to bother with adjusting illumination. The display, with its ambient light sensor, automatically adjusts screen brightness to maintain readability, in nearly any lighting environment.

The monitor also has built-in Advanced Clinical Solutions that provide tools to summarize and visualize complex clinical data and their interactions. Multiple streams of information come together in one, uncomplicated view for you.



Ambient light adjustment is just part of why the IntelliVue MX550 is an excellent solution for critical care settings and the OR. In the OR, for example, the display automatically adapts to lights being dimmed or brightened during different stages of surgery, reducing overall light pollution in the OR.



Make the most of your EMR

EMRs represent a considerable financial expense. The IntelliVue MX550 works to enhance your investment by sharing data with clinical information systems helping to contribute to a more comprehensive EMR. It can also help reduce cost and complexity when connecting bedside devices to your EMR solution of choice. Philips IntelliBridge device interfacing makes this possible by eliminating the need for a separate device concentrator and data consolidation server.

IntelliBridge technology* provides a single point of interfacing for ventilators and other devices at the bedside and during in-hospital transport. You can add up to two device interfaces as expansion boards, IntelliBridge EC10 modules, or a combination of the two.

Clinical data from the monitor and other linked bedside devices are combined into a single HL7 message stream** for reliable, standards-based interfacing to EMRs. Even alarm strips captured as electronic documents during transport can be transferred using a document import interface provided by the EMR. The IntelliVue MX550 can thus reduce the need for expensive bedside device concentrators and associated licensing contracts.

* Requires IntelliBridge interface port(s). For device compatibility, see the latest IntelliBridge EC10 external device compatibility list.

** Wireless data transfer to EMR via Philips IntelliVue Information Center (PIIC) or Philips IntelliVue Information Center iX (PIIC iX)



Fit monitoring around workflows

Today's workflows need to adapt in the blink of an eye. The IntelliVue MX550 can be scaled to cover most critical care monitoring requirements as levels of patient acuity change. It supports the Philips Multi Measurement Module (MMS), IntelliVue X2 module, Measurement Server Extensions (MMSE), and specialty measurements via up to three single/double-width parameter modules as required.



This monitor is designed to simplify clinical workflows. Since the IntelliVue MX550 can interface to other bedside devices, there's no need for a bedside device concentrator or an interface server.

Take monitoring with you

The IntelliVue MX550 delivers the functionality of a full-fledged bedside monitor while providing battery back-up to bridge the occasional power loss. With standard battery operation, a built-in handle, and overall rugged design, this solution can cope with the demanding in-hospital transport sometimes required in acute care areas. It's also compact enough to be convenient underway. If transport is necessary, the optional Philips quick-release mount allows you to unfasten this monitor from a secure bedside mounting in seconds. You can also keep the pace moving with the optional and lightweight bed hanger mount developed exclusively for the IntelliVue MX series.

Lower the likelihood of cross-contamination

We've engineered the IntelliVue MX550 with infection control in mind. It's simple to clean thanks to smooth surfaces, minimal seams, and sturdy material approved for use with a variety of tested disinfecting agents.¹

As with all IntelliVue monitors, this monitor has no hard drive and no cooling fans. As fan blades are likely to accumulate and distribute dust and pathogens, the IntelliVue MX550 features passive cooling, helping you in your battle against nosocomial infections.

¹ Refer to the IntelliVue Patient Monitoring Instructions for Use for a detailed list of approved cleaning agents.



**Philips Healthcare is part of
Royal Philips**

How to reach us

www.philips.com/healthcare

healthcare@philips.com

Please visit www.philips.com/IntelliVueMX550



© 2014 Koninklijke Philips N.V.
All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Printed in The Netherlands
4522 991 01811 * FEB 2014