
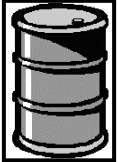




Philips Medical Systems/Product Recycling Passport

Printed copies are uncontrolled

Product name:	M3160A 4-channel Thermal Recorder	
Identification code(s)	M3160-60001, M3160A	
Total weight (in Kg)	2.9	
Producer/ Manufacturer	Name company:	Philips Medical Systems
	Address:	3000 MinuteMan Road, Andover, MA
	Zip code:	01810
	Country:	USA

Recycle Info	Items:	Location
Special attention 	Remove Bottom screws on unit and open front panel	Figure 1 & 2
Fluids / Gases 	Items:	Location
	Not Applicable	
Batteries 	Type:	Location
	Not Applicable	
To be Removed		
Hazardous 	Substances:	Location
	2 Printed Circuit Boards	Figure 1 & 3
To be Removed		

Title: M3160A 4-Channel Thermal Recorder Recycling Passport
 Recycling Passport Number: A-M3160-90002 Rev: A

Philips Medical Systems/Product Recycling Passport

Printed copies are uncontrolled

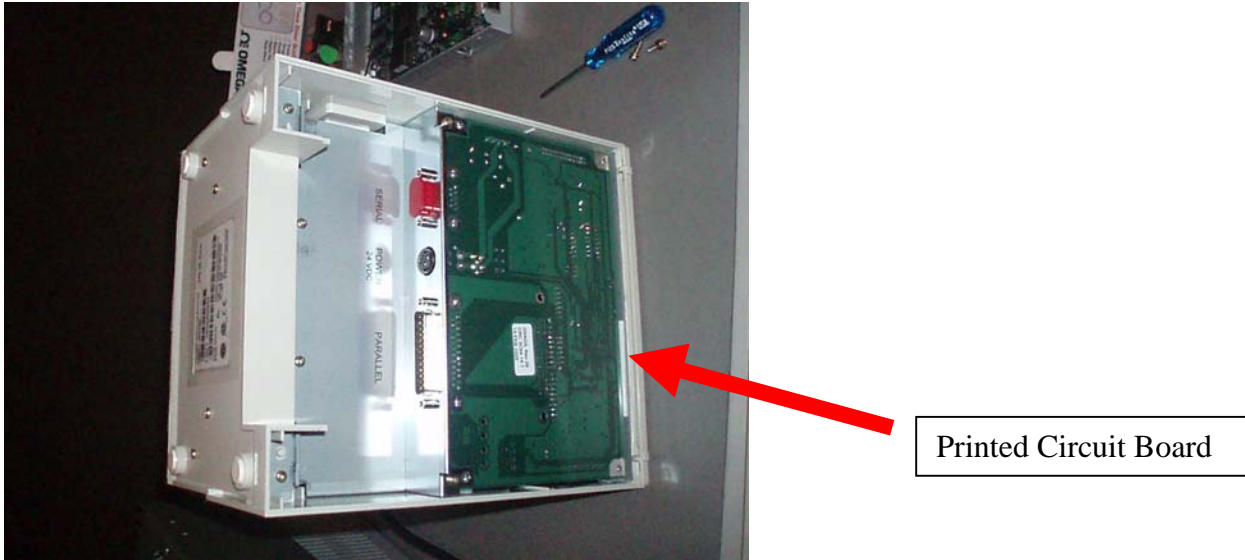


Figure 1 (Remove screws on bottom of unit with Phillips Screwdriver)



Figure 2 (Open front panel)

Title: M3160A 4-Channel Thermal Recorder Recycling Passport
Recycling Passport Number: A-M3160-90002 Rev: A

Philips Medical Systems/Product Recycling Passport

Printed copies are uncontrolled

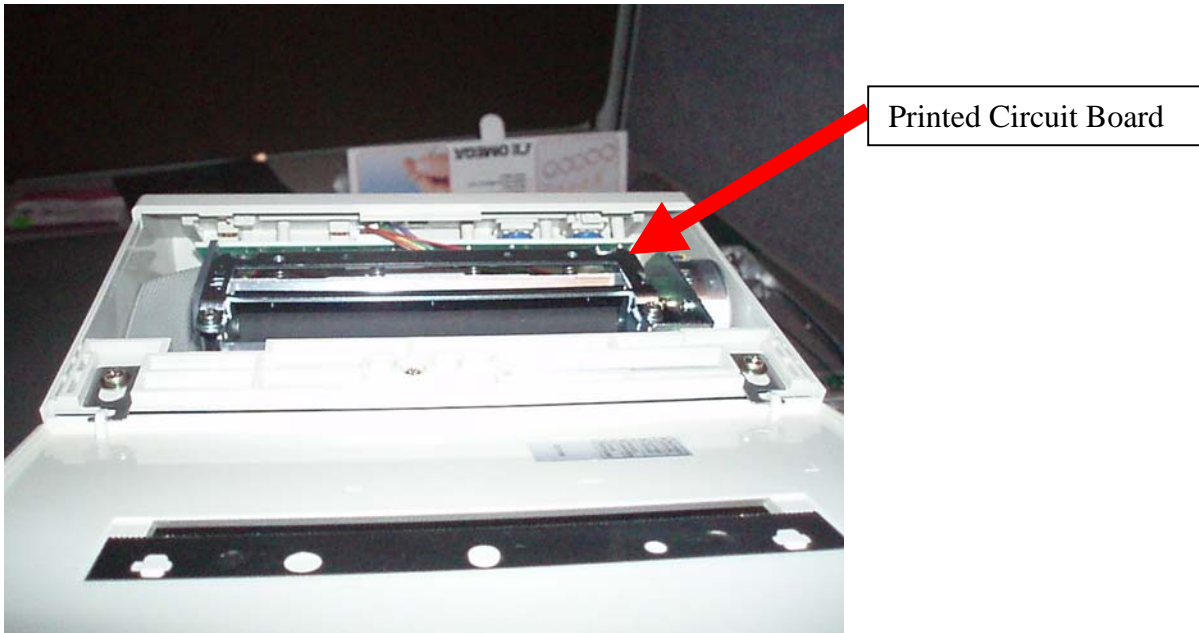


Figure 3 (rotate unit and look up under LEDs and touch buttons)