
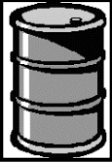




<b>Product name:</b>	<b>ClearVue 550 Medical Diagnostic Ultrasound Equipment</b>	
<b>Identification code(s)</b>	<b>989605409201</b>	
<b>Total weight (in Kg)</b>	<b>56.43 kg</b>	
<b>Producer/ Manufacturer</b>	<b>Name company:</b>	Philips Ultrasound
	<b>Address:</b>	22100 Bothell Everett Highway, Bothell, WA
	<b>Zip code:</b>	98012
	<b>Country:</b>	USA

Recycle Info	Items:	Location
<b>Special attention</b>  	Gas strut under pressure	Figure 1 [1]

Fluids / Gases	Items:	Location
	Oil and Nitrogen gas: in gas strut in Articulation Arm mechanism. Do not drain oil with device under pressure	Figure 1 [1]

Batteries	Type:	Location
  <b>To be Removed</b>	Battery in System. Button cell on Motherboard.	Figure 1 [2]

Hazardous	Substances:	Location
  <b>To be Removed</b>	Mercury (Hg) – Monitor contains mercury in the backlight lamps.	Figure 1 [3]

Show locations of materials mentioned on the previous sheet. (Picture information)

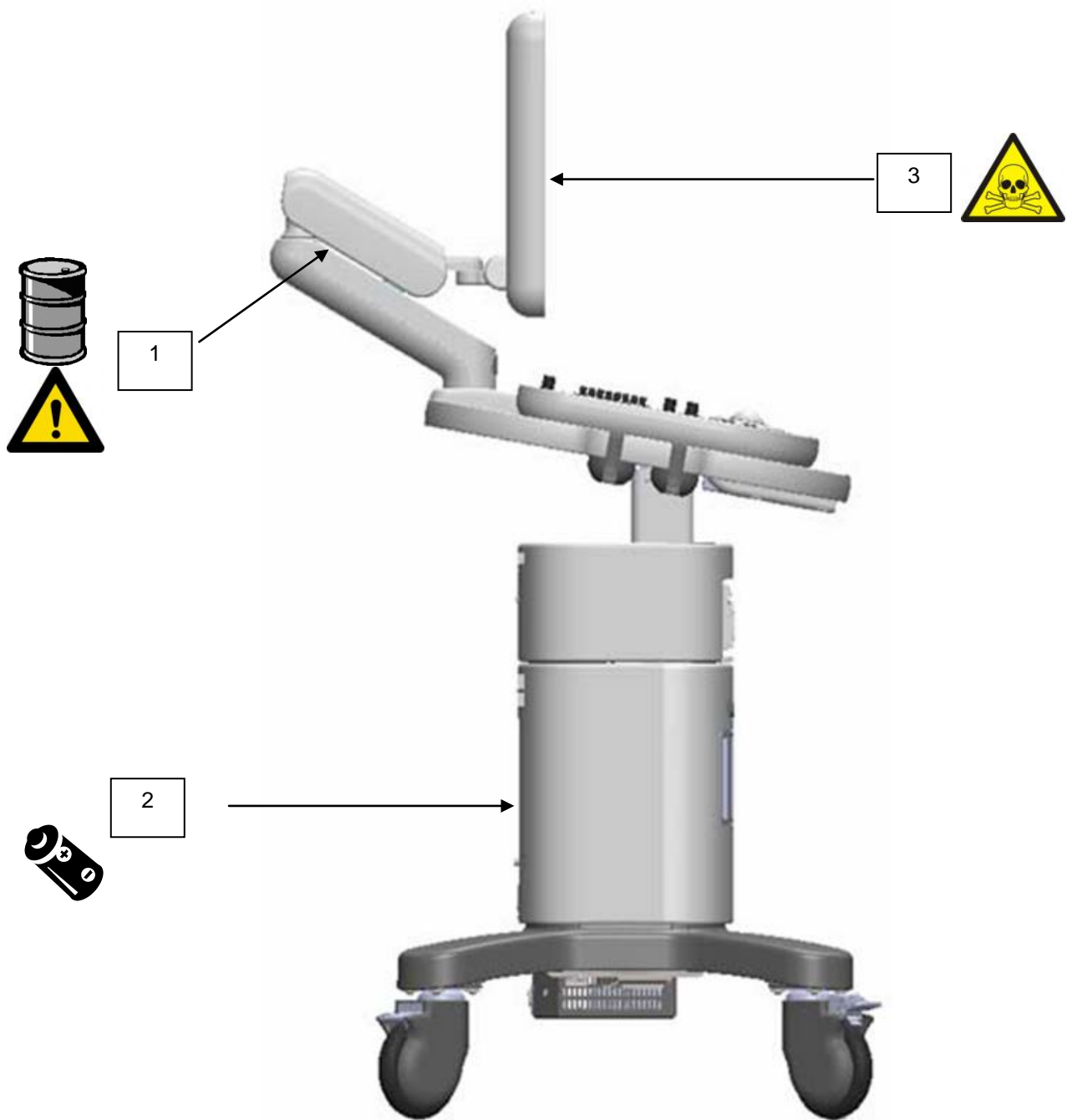


Figure 1: ClearVue 550