

<b>Product name:</b>	<b>EPIQ7 Medical Diagnostic Ultrasound Equipment</b>	
<b>Identification code(s)</b>	<b>989605386721</b>	
<b>Total weight (in Kg)</b>	<b>105 Kg</b>	
<b>Producer/ Manufacturer</b>	Name company:	Philips Ultrasound
	Address:	22100 Bothell Everett Highway, Bothell, WA
	Zip code:	98021
	Country:	USA


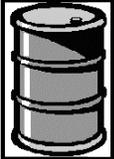
Recycle Info	Items:	Location
<b>Special attention</b> 	Gas Spring under pressure.	Figure 1 [1], [2]
<b>Fluids / Gases</b> 	<b>Items:</b>	<b>Location</b>
	Oil and Nitrogen Gas: (2x) in Gas Spring Tubes in Control Panel Articulation Assy Do not drain oil with device under pressure.	Figure 1 [1]
	Oil and Nitrogen Gas: in Gas Spring Tube in Arm Articulation Assy Do not drain oil with device under pressure.	Figure 1 [2]
<b>Batteries</b>  <b>To be Re moved</b>	<b>Type:</b>	<b>Location</b>
	Li Coin cell (1x) on CPU PCB	Figure 1 [3]
	Batteries in System (2x) Li-Ion batteries in power system	Figure 1 [4]
<b>Hazardous</b>  <b>To be Removed</b>	<b>Substances:</b>	<b>Location</b>
	Lead (Pb) - contains lead in PCB's, solder, wiring, and cables	N/A
	Mercury (Hg) – Video monitor lamps	N/A

Figure 1

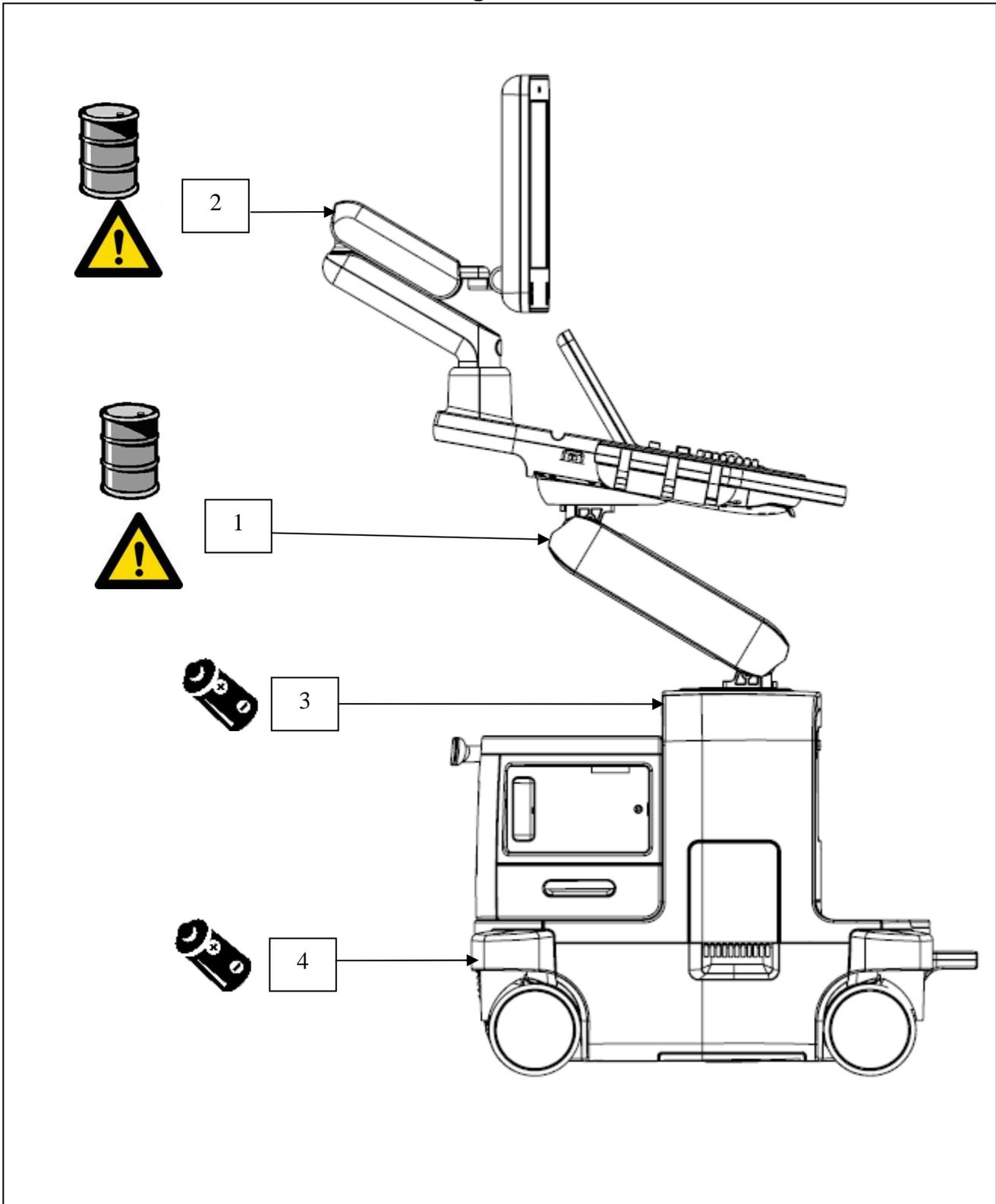

**Recycling Passport Revision History**

<b>Rev</b>	<b>Description</b>	<b>Date</b>
B	Revision for Released Product Name, EPIQ7	July 24, 2013

**\*\*\*\*\*End of Document\*\*\*\*\***